**NATURE &
DECOUVERTES**

MASSEUR RÉFLEXOLOGIE
REFLEXOLOGY MASSAGER
MASAJEADOR DE REFLEXOLOGIA
MASSAGEAPPARAAT VOOR VOETREFLEXOLOGIE
MASSAJADOR DE REFLEXOLOGIA
Réf. 15229380

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PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

IMPORTANT SAFEGUARDS

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning. To disconnect, turn all controls to the 'OFF' position, then remove plug from outlet.
- NEVER leave an appliance unattended when plugged in. Unplug from the mains outlet when not in use and before putting on or taking off parts or attachments.
- This appliance is designed for use by adults. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by unsupervised children.
- DO NOT reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug it immediately. Keep dry – DO NOT operate in wet or moist conditions.
- NEVER insert pins or other metallic fasteners into the appliance.
- Use this appliance for the intended use as described in this booklet.
- DO NOT use attachments not recommended by the manufacturer.
- NEVER operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the Point of Purchase for examination and repair.
- Keep cord away from heated surfaces.
- Excessive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by the cord or use cord as a handle.
- This appliance is for household and indoor use only.
- DO NOT crush.
- This product requires a 220-240 V AC power.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. For service, return to the Point of Purchase. Servicing of this appliance must be performed by authorised service personnel only.
- Do not stand on or in the device. Use only while seated.
- Never block the air openings of the appliance or place it on a soft surface, such as a bed or couch. Keep air openings free of lint, hair etc.

- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- Failure to follow the above may result in the risk of fire or injury.
- Burns can result from improper use.

SAFETY CAUTIONS


PLEASE READ THIS SECTION CAREFULLY BEFORE USING THE APPLIANCE.

- If you have any concerns regarding your health, consult a doctor before using this appliance.
- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance. Not recommended for use by diabetics.
- DO NOT use on an infant, persons with physical frailty, or on a sleeping or unconscious person. DO NOT use on sensitive skin or on a person with poor blood circulation.
- This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the lower half of their body.
- NEVER use directly on swollen or inflamed areas or skin eruptions.
- If you feel any discomfort whilst using this appliance, discontinue use and consult your doctor.
- This is a non-professional appliance, designed for personal use and intended to provide a soothing massage to worn muscles. DO NOT use as a substitute for medical attention.
- DO NOT use before going to bed. The massage has a stimulating effect and can delay sleep.
- Do not use for longer than the recommended time.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- Care should be taken when using heated surfaces. If the product feels excessively hot, switch off at mains and contact the Point of Purchase.
- Please ensure that all hair and jewellery is kept clear of the massage mechanism, or any other moving parts of the product, at all times.
- Take care when using and storing the appliance and its power cable so that it is not a trip hazard.
- To prevent trips, take care to fully remove both feet from the device before movement or activity.
- The device must not be used by infants (0-3 years).

STORAGE

Place the appliance in its box or in a safe, dry, cool place. To avoid breakage, DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.

PRODUCT OVERVIEW

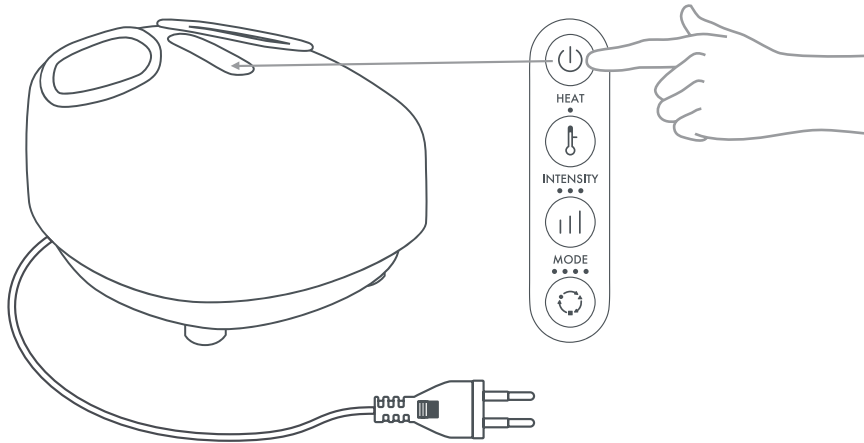POWER CONNECTION / MAINS PLUG

Power button
Turns the device on/off

Heat
Toggles the warming mode on/off

Intensity Level
(low, medium, high)


Mode select (1,2,3,4)


SET-UP

Before using the massager:

- Connect it to a 220-240v mains outlet
- Place on the floor in front of a chair
- We recommend that you wear socks for optimum comfort
- Sit down, then place your feet into the two openings on the device.

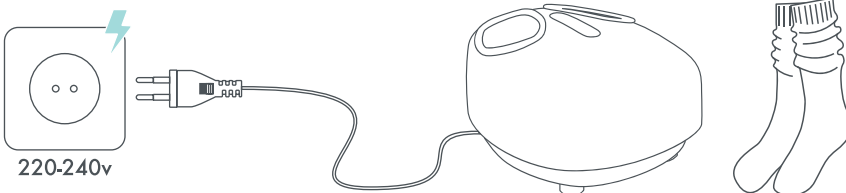
TURNING THE POWER ON

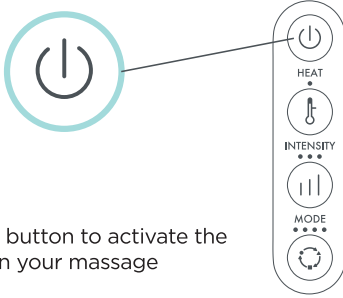Press the power button to activate the device and begin your massage session.


HEAT

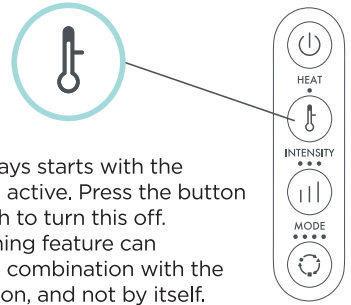The device always starts with the warming mode active. Press the button once if you wish to turn this off. Note: the warming feature can only be used in combination with the massage function, and not by itself.

AIR MASSAGE - INTENSITY


Press the intensity button to cycle through the 3 intensity levels (low, medium, high) to find the one that is most comfortable and effective for you. The LED lights indicate which level you currently have selected. The device always starts at intensity level 1, the lowest level of intensity.

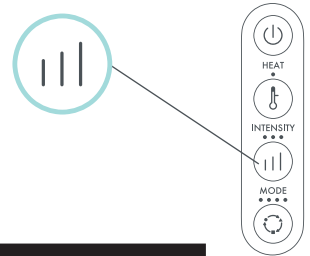
MODE SELECT

While the device is active, press this button to cycle through 4 pre-programmed massage modes. The LED lights indicate which mode you currently have selected.


Mode 1 : focused air pressure massage.
With secondary rolling massage.

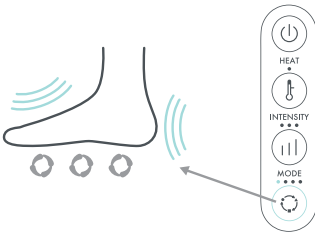
Mode 2 : focused rolling massage.
With secondary air pressure massage.

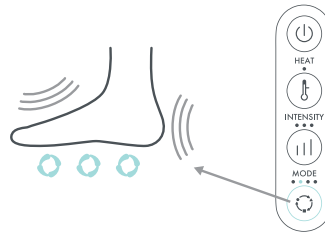Mode 3 : air pressure massage only

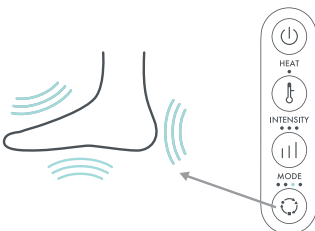Mode 4 : Rolling massage only

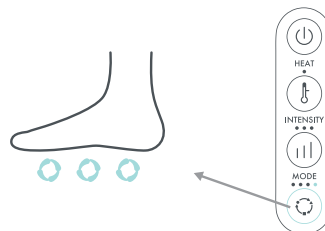SHUTDOWN

Each session lasts for 15 minutes. The device automatically shuts down after this time. At any other time during your massage session you can press the power button to turn off the device as required.


15
minutes


To Restart

WARNINGS AND GUIDANCE


- ⚠️ Only use the massager while seated. Do not stand in or on the massager.
- Always ensure the device is turned off when not in use and unplugged from the mains power
- Do not use the device continuously for longer than 15 minutes. Take a break between sessions.

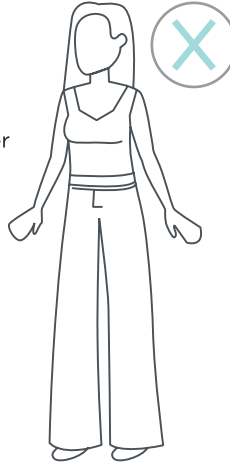MAINTENANCE

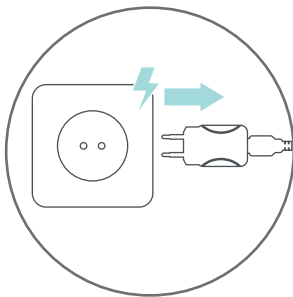Disconnect from power source before cleaning.

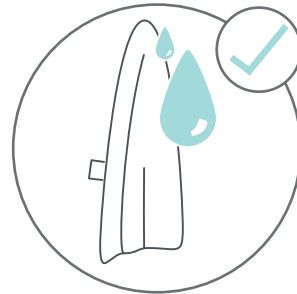May be wiped with a soft, very lightly damp cloth or sponge.

Avoid bringing the device into contact with liquids and never immerse the device in liquid.

Do not apply cleaning products to the device, or any abrasive cleaning appliances

SPECIFICATION

- Model: reflexology massager - ZEN-F1-MP
- Product Dimensions: L 445 mm x W 408 mm x H 231 mm
- Rated Voltage: 220V-240V-, 50-60Hz
- Rated Power: 50W
- Time Rating: 15 minutes
- Noise: ≤ 60 dB

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