
**MANUEL D'UTILISATION
USER MANUAL
BENUTZER-HANDBUCH**

**MASQUE DE RELAXATION
EYE MASSAGER
MASK ENTSPANNUNG**

Réf. 15193130

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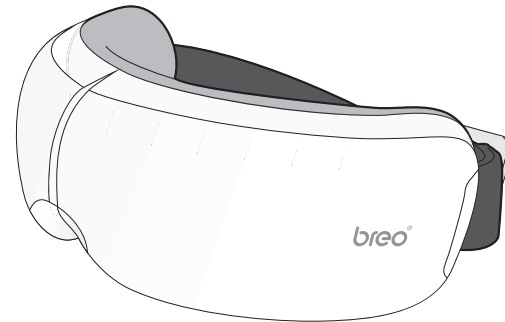11. It is household electric health care product. thus can not be used for diagnosing or curing known disease or other commercial use.
12. If it requires cleaning, wipe it with cotton soaked in alcohol. Do not use water to wash the device.
13. The massager has a hot compressing function, and temperature can reach at 42°C, users not sensitive to change in temperature should be careful in use.
14. Take off your eye contact lenses or glasses before use this product.

TROUBLESHOOTING AND TECHNICAL PERFORMANCE INDEX

If there is any problem when you use the device, please follow the instructions below to check. If the problem still exists or the device can't work normally, please turn off the power supply and contact local service center. Do not disassemble or repair the device yourself.

Fault	Possible reason	Solutions
Unable to boot	Whether the cable is in good condition	Check the connection status of the cable
	Whether the battery is low	If yes, charge battery
Air pressure is too weak.	Whether the battery is low	If yes, Please charge at once

TECHNICAL PERFORMANCE INDEX

Input	DC:5V 1A
Power	5W
Ideal environment	Temperature: 5°C - 40°C - Relative Humidity: <80% Atmospheric pressure: 860hpa-1060hpa
Dimension	Master unit: L.208 mm x W.106 mm x H.76 mm
Net Weight	318 g

OPERATION

Mode Selection

After starting up, the machine's default mode is SLEEP mode; short-press On/Off key to switch next mode.

SLEEP MODE

Time: 10 minutes
Airpressure: ON
Hot compressing: ON
Vibration: OFF
Music: OFF

MED MODE

Time: 15 minutes
Airpressure: ON
Hot compressing: ON
Vibration: OFF
Music: ON

HARD MODE

Time: 20 minutes
Airpressure: ON
Hot compressing: ON
Vibration: ON
Music: ON

Shutdown

Long-press On/Off key to shut down the machine, and the machine stops working.

PRECAUTIONS

Please read User Manual carefully to correctly and safely use this product.

1. Please do not use this product while driving.
2. Please do not use this product in the environment with excessively high temperature or humidity.
3. Wear the device and then switch it on to avoid damaging the air bag due to overpressure.
4. Those suffering from eye diseases. such as glaucoma. pearl eye and detachment of retina can not use this product
5. Those who have had eye surgery can not use this product.
6. The persons suffering from strokes. hypertension. and hypotension should not use this product.
7. Sluggish patients, or people with physical or mental or intellectual disability should not allowed to use this product
8. Those who have eye surgery can not use this product.
9. Children under 7 years old cannot use this product.
10. Please use the stabilized voltage power supply provided by our company and pay attention to the safety.

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
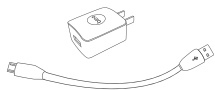
Handbuch

FUNCTION CONFIGURATION


1. Power Supply Indicator
2. On/Off Key
3. Display Area
4. Elastic Adjustment Buckle
5. USB Input Socket

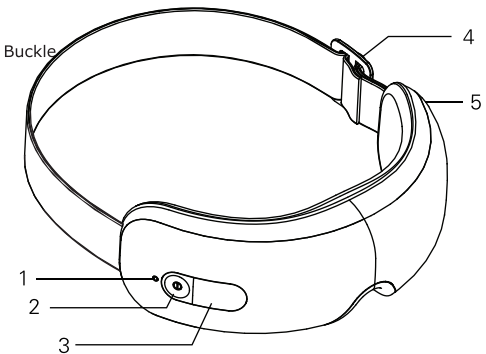OPERATION

Power on

1. Open the lid of the charging port in the bottom of the rear cover of the master unit.
2. Plug the power cord into a wall outlet. We suggest using a surge protected power bar, and plug the other end into the USB jack under the rear cover.
3. The indicator of power on the LCD screen flashes during charging; the flash stops when fully charged.
4. The product cannot be operated while charging (if the product is charged while operating, it will start the self-protection program and stop operating automatically)
5. The time of charging is 2 -3hours.

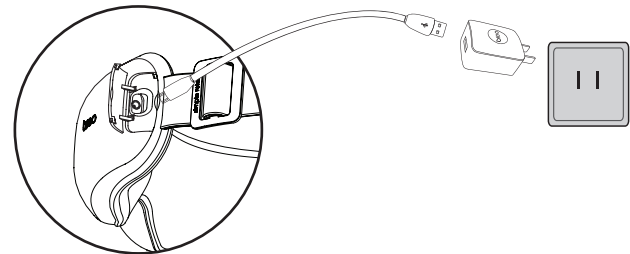USER MANUAL

EYE MASSAGER

Ref. 15193130


Thank you for purchasing Breo iSee4
Please read these instructions carefully before using and keep it properly for future reference.

INTRODUCTION

iSee4 eye massager adopts intelligent air pressure, micro-computer chip control technology. By kneading, hot compressing and massaging on eyes, it helps to keep your eyes in the healthiest state, and while performing massage with this product, you can enjoy the wonderful music. With the unique 180°-foldable and integrated design, the product can fit eyes completely, and make people full of vim and vigor.

PRODUCT CONFIGURATION

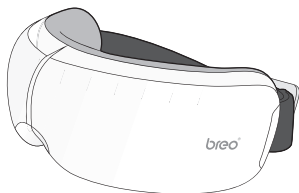The massager includes : Master Unit, Power Charger, User Manual.

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Copyright © Breo. All Rights Reserved.