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### **1. About the SAD Light**

#### **Instructions for Use**

Congratulations on the purchase of your new “SAD Light”. The product will offer many years of reliable usage if handled correctly. Please keep your receipt, these instructions and the original packaging safe in case of future requirement.

Please read the instructions carefully to gain the most benefit from your new “SAD Light”, although the usage guidelines in this document do not replace the advice of your qualified medical advisor.

#### **Intended Use:**

Medical intended use of this appliance is to treat the symptoms SAD (seasonal affective disorder) or winter depression. “SAD Light” is intended to boost people energy level, to adjust circadian rhythms and make people feel more awake. “SAD Light” is safe and medically proven to provide relief from the symptoms of SAD.

#### **How does the “SAD Light” work?**

In autumn and winter, the seasons with the least sunlight because of shorter days, the lack of sunlight and a problem with certain brain chemicals stops the hypothalamus working properly. The lack of light is thought to affect the production of the hormone melatonin (that makes you sleepy).

“SAD Light” uses bright lights to simulate the sunlight during the darker autumn and winter, whilst eliminating the dangerous UV rays contained in sunlight. The additional light helps to adjust the levels of melatonin that can help the symptoms of SAD.

#### **Seasonal Affective Disorder:**

Seasonal Affective Disorder (SAD) as known as winter depression or winter blues. It may cause people:

- Feel low, listless and fatigued
- Have trouble concentrating
- Lack in energy
- Need for more sleep
- Less active than normal

Symptoms occur during the same season each year. Most people with SAD start to have symptom in winter and feel better in summer.

#### **Contraindication:**

The following conditions may be contraindication for use of this appliance:

- suffer(ed) from severe depression
- recent eye surgery or a diagnosed eye condition whereby doctor advised you to avoid bright light
- on specific medication (i.e. certain antidepressants, psychotropic drugs or malaria tablets);
- have high blood pressure with vitreous hemorrhage have no benefit;

- skin sensitive to light, such as systemic lupus erythematosus

**Always consult your doctor before your start using this appliance if you have any of the above conditions.**

**Possible side effects:**

The following side effects you may experience:


- headaches
- Eyestrain
- sleeping problems
- Nausea

Measures to avoid, minimize or alleviate these side effects:

- Increase the treatment distance
- Avoid using the appliance during the late evening or before sleep
- Stop using the appliance for a few days to make the side effects disappear and then try again.

## 2. Signs and symbols

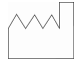The following symbols appear in this manual, on the “SAD Light” or on its accessories:

 <b>Warning</b>	Warning instruction indicating a risk of injury or damage to health		Read operating instructions before use
 <b>Important</b>	Safety note indicating possible damage to the appliance/ accessory		Caution: Consult accompanying documents
<b>SN</b>	Specifies serial number of the “SAD Light”		Class II Equipment
	OFF/ ON		Date of Manufacture
	Authorized Representative in the European Community		Manufacturer
<b>IP20</b>	Ingress Protection Rating		Keep Dry

## 3. Safety and performance standards

The “SAD Light” has been designed and manufactured to conform to the highest standards of safety and performance including electromagnetic compatibility (EMC). The “SAD Light” conforms to the applicable requirements of the following:

This appliance meets the requirement of Medical Device Directive 93/42/EEC for medical appliances.

## 4. Notes

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

### Safety notes

#### **Warning**

- Do not use this appliance in wet surroundings (e.g. in the bathroom or near a shower). Do not immerse appliance in water or any other liquid.
- Before use, ensure that all packaging materials are removed and that there is no visible damage to the appliance or accessories. When in doubt, do not use the appliance and contact your dealer or the customer service address provided.
- The adaptor contains a transformer. Do not cut off the adaptor to replace it with another plug, as this can cause a hazardous situation.
- This appliance is only intended for use on the human body.
- Close supervision is necessary when this appliance is used by, or near children or invalids.
- Children should be supervised to ensure that they do not play with the appliance.
- There is no need for any special training for this appliance.
- Keep the cord away from hot objects and open flames.
- Connect the appliance only to the mains voltage listed on the type plate.
- Always remove the mains plug and let the appliance cool off before you touch it.
- Do not touch the appliance with wet hands while it is plugged in; do not allow any water to be sprayed onto it. The appliance must be operated only when it is completely dry.
- Protect the appliance from strong impacts.
- Avoid pulling the cord when removing the plug from the socket. Instead, hold the plug.
- If the adaptor, cord or the appliance is damaged, contact customer service or the dealer, since special tools are needed to make repairs. Never try to repair the appliance yourself.
- Disconnection from the power supply network is only guaranteed when the mains part is unplugged.
- This appliance is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.
- This appliance requires no calibration, preventive inspection and maintenance.
- This appliance is not repairable and contains no user serviceable parts.
- The operator shall not touch the appliance and the patient simultaneously.
- The appliance is protected against unauthorized use. Do not modify this appliance without authorization of the manufacturer.
- If this appliance is modified, appropriate inspection and testing must be conducted to ensure continued safe use of the equipment.

### General notes

#### **Important**

- This appliance is not intended for commercial or clinical use, but only for individual, private household use.
- If you are taking medication such as pain relievers, antihypertensives or antidepressants, consult your doctor before using this appliance.
- People with retinal diseases as well as diabetics need to be examined by an ophthalmologist before using this appliance.
- Please do not use if you have an eye disease such as cataracts, glaucoma, diseases of the optic nerve in general and inflammation of the vitreous humour.
- If you have health concerns of any kind, contact your general practitioner!

## 4. Notes (Continued)

### EMC information

#### **⚠ Warning**


- The appliance needs special precaution regarding EMC and needs to be installed and put into service according to the EMC information provided in the section “Normative reference”.
- The appliance needs to be installed and put into service in accordance with the information provided in the accompanying documents.
- Wireless communications equipment such as wireless home network devices, mobile phones, cordless telephones and their base stations and walkie-talkies can affect this appliance and should be kept at least a distance (d) away from the appliance. The distance (d) is calculated by the manufacturer from the 800 MHz to 2.5 GHz column of table in the section “Normative reference.

### Repairs

#### **⚠ Important**

- Please do not attempt to repair the appliance yourself. This could result in serious injuries. Failure to observe this regulation shall void the warranty.
- If this appliance is damaged in any way or it is not working, it should not be operated. Please try to fix the problem according to FAQ table. If the problem is still not resolved, please contact customer service or the dealer, since special tools are needed to make repairs. Never try to repair the appliance yourself.

## 5. Unit description

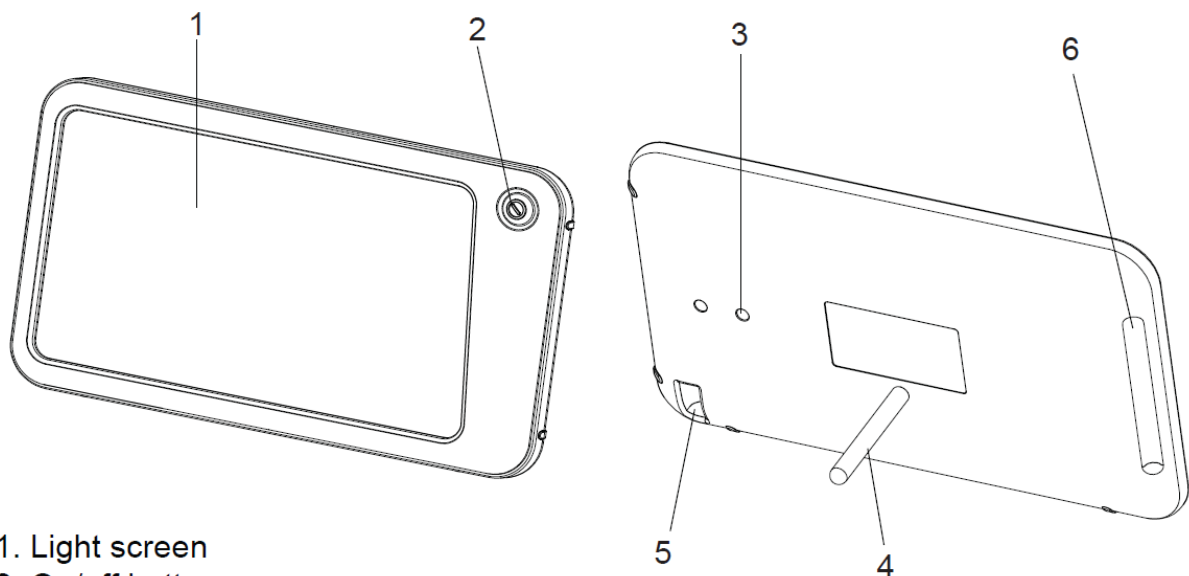1. Light screen
2. On/off button
3. Stand hole
4. Stand
5. Socket for appliance plug
6. Stand storage cavity

## 6. Installation

1. Remove the stand from the stand storage cavity.
2. Insert the end of stand into the hole on the back of appliance.
3. Place the appliance on a solid surface approximately 20cm from your eyes, the light must bathe your face from the side.
4. Connect the appliance only to the mains voltage listed on the type plate.

### Note

- Push the plug completely into the mains socket.
- Ensure that there is a socket near the place where the appliance will stand.
- Lay the mains cable in such a way that no one can trip over it.

## 7. Operation

### 1. Turning on

Press the ON/OFF button once of the appliance to turn on.

### 2. Taking the light

You can read, write, make telephone calls or work while you use the appliance. When using the appliance, try to maintain the recommended distance of 20 cm between your eyes and the appliance. The duration of the treatment is approx. 2 hours.

### Note

- Do not look directly into the light continuously! Some eye pain or mild headaches may be experienced during initial sessions but these should pass. If in doubt just lower your treatment times until comfortable again.
- Do not leave the appliance unattended in the duration.

### 3. Enjoying light over longer periods

During the darkest period of the year, repeat the treatment for at least 7 consecutive days, or longer, depending on your individual requirements.

### 4. Switching off

Press the ON/OFF button of the appliance to switch off.

## 8. Cleaning

### Warning

- Before you clean the appliance, always switch it off, unplug it and let it cool off.
- Do not touch the appliance with wet hands.
- Clean using a damp cloth, to which you can apply a little detergent if necessary.
- Do not let water run into the appliance or spill water onto the appliance.
- Cleaning can be done once per month.

## 9. Storage

When this appliance is not being used, it should be stored in a dry location out of reach of children. Do not keep the appliance connected to a power source. Follow the storage instructions provided in the section “Technical Specification”.

## 10. Disposal

Please note that this appliance is marked with the WEEE symbol on the left. This means that this appliance must not be disposed of together with ordinary

household waste, as electrical and electronic waste must be disposed of separately.

If you have any questions, please contact your local authorities responsible for waste disposal.

## 11. FAQ

Problem	Possible causes	Solutions
Device does not light up	On/ Off button switched off	Switch on the on/off button
	No Power	Connect the main plug correctly
	No Power	The main cable is defective. Contact customer services or your retailer
	LEDs have reached the end of its service life	If you need to have the appliance repaired, contact customer service or an authorized dealer.
	LEDs faulty	

## 12. Technical specification

<b>Model:</b>	LND018915V3
<b>Dimension (WHD)</b>	260 x 150 x 17 mm
<b>Weight</b>	450g
<b>Lighting elements</b>	LEDs
<b>Power</b>	DC12V, 1A
<b>Power</b>	12 watts
<b>Light intensity</b>	> 10,000 lux in 6 cm
<b>Mains connection</b>	12V DC, 1.0A, only in adaptor model CS12F120100FGF or CS12F120100FBF
<b>Operating Environment</b>	0 to +40 °C, and 30% to 90% relative humidity Atmospheric pressure: 700-1060hPa
<b>Storage Environment</b>	-10 to +60 °C, and 30% to 90% relative humidity Atmospheric pressure: 700-1060hPa
<b>Transport Environment</b>	-10 to +60 °C, and 30% to 90% relative humidity Atmospheric pressure: 700-1060hPa
<b>Product Classification</b>	1. Class II equipment 2. No applied part 3. Protection against ingress of water: IP20 4. Not category AP / APG equipment 5. Mode of operation: Continuous operation.
<b>Maximum output of radiation/ spectral radiant</b>	14 W.m <sup>-2</sup> .sr <sup>-1</sup>
<b>Location of emission aperture</b>	Screen
<b>Treatment area</b>	Around 400 cm <sup>2</sup>
<b>Risk Group of light source equipment</b>	Exempt group classified by IEC 60601-2-57:2011
<b>Product Service Life</b>	5 years

Subject to technical modification.

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### **European Authorised Representative:**

I-Plus Company Ltd.  
4B Christchurch House, Beaufort Court,  
Medway City Estate, Rochester, Kent. ME2 4FZ, UK  
Tel: (44)-20-71936086  
Email: i-plus@zeroplus.hk

### **Important**

If this appliance is damaged in any way or it is not working and you cannot find the answer to your question in FAQ. Please contact customer service or the retailer for assistance.

## 14. Normative references

<b>Guidance and manufacturer's declaration – electromagnetic immunity</b>			
The model LND018915V3 is intended for use in the electromagnetic environment specified below. The customer or the user of the model LND018915V3 should ensure that it is used in such an environment.			
<b>Immunity test</b>	<b>IEC 60601 test level</b>	<b>Compliance level</b>	<b>Electromagnetic environment – guidance</b>
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines	±2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV differential mode	±1 kV differential mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5 % $U_T$ (>95 % dip in $U_T$ ) for 0,5 cycle  40 % $U_T$ (60 % dip in $U_T$ ) for 5 cycles  70 % $U_T$ (30 % dip in $U_T$ ) for 25 cycles  <5 % $U_T$ (>95 % dip in $U_T$ ) for 5 sec	<5 % $U_T$ (>95 % dip in $U_T$ ) for 0,5 cycle  40 % $U_T$ (60 % dip in $U_T$ ) for 5 cycles  70 % $U_T$ (30 % dip in $U_T$ ) for 25 cycles  <5 % $U_T$ (>95 % dip in $U_T$ ) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the model LND018915V3 requires continued operation during power mains interruptions, it is recommended that the model LND018915V3 be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE $U_T$ is the a.c. mains voltage prior to application of the test level.			

## 14. Normative references (Continued)

<b>Guidance and manufacturer's declaration – electromagnetic immunity</b>			
The model LND018915V3 is intended for use in the electromagnetic environment specified below. The customer or the user of the model LND018915V3 should ensure that it is used in such an environment.			
<b>Immunity test</b>	<b>IEC 60601 test level</b>	<b>Compliance level</b>	<b>Electromagnetic environment – guidance</b>
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 V	<p>Portable and mobile RF communications equipment should be used no closer to any part of the model LND018915V3, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p><b>Recommended separation distance</b></p> <p><b><math>d=1.167\sqrt{P}</math></b></p> <p><b><math>d=1.167\sqrt{P}</math> 80MHz to 800MHz</b></p> <p><b><math>d=2.333 \sqrt{P}</math> 800MHz to 2.5 GHz</b></p>
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2,5 GHz	3 V/m	<p>where <math>P</math> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and <math>d</math> is the recommended separation distance in metres (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,<sup>a</sup> should be less than the compliance level in each frequency range.<sup>b</sup> Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
<p>NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.</p> <p>NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p>			
<p><sup>a</sup> Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the model LND018915V3 is used exceeds the applicable RF compliance level above, the model LND018915V3 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the model LND018915V3.</p> <p><sup>b</sup> Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [V1] V/m.</p>			

## 14. Normative references (Continued)

<b>Guidance and manufacturer's declaration – electromagnetic emissions</b>		
The model LND018915V3 is intended for use in the electromagnetic environment specified below. The customer or the user of the model LND018915V3 should ensure that it is used in such an environment.		
<b>Emissions test</b>	<b>Compliance</b>	<b>Electromagnetic environment – guidance</b>
RF emissions CISPR 11	Group 1	The model LND018915V3 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The model LND018915V3 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class C	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

<b>Recommended separation distances between portable and mobile RF communications equipment and the LND018915V3</b>			
The model LND018915V3 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Model LND018915V3 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the model LND018915V3 as recommended below, according to the maximum output power of the communications equipment.			
<b>Rated maximum output power of transmitter W</b>	<b>Separation distance according to frequency of transmitter m</b>		
	<b>150 kHz to 80 MHz</b> $d = \left[\frac{3,5}{V_1}\right]\sqrt{P}$	<b>80 MHz to 800 MHz</b> $d = \left[\frac{3,5}{E_1}\right]\sqrt{P}$	<b>800 MHz to 2,5 GHz</b> $d = \left[\frac{7}{E_1}\right]\sqrt{P}$
0,01	<b>0.117</b>	<b>0.117</b>	<b>0.233</b>
0,1	<b>0.369</b>	<b>0.369</b>	<b>0.738</b>
1	<b>1.167</b>	<b>1.167</b>	<b>2.333</b>
10	<b>3.690</b>	<b>3.690</b>	<b>7.377</b>
100	<b>11.67</b>	<b>11.67</b>	<b>23.33</b>
For transmitters rated at a maximum output power not listed above, the recommended separation distance $d$ in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where $P$ is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

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