NATURE &  
DECOUVERTES

# MASSEUR CERVICAL BREO

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## NECK MASSAGER BREO

# PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE.

## WARNING


Please carefully read this manual and keep it in a safe place

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- The appliance has a heated surface.
- Persons insensitive to heat must be careful when using the appliance.
- Do not insert any object or flammable material into the neck and shoulder massager. This may result in electric shock, malfunction or risk of fire.
- Do not exert excessive pressure when using the neck and shoulder massager.
- Do not cover the neck and shoulder massager during operation. This may result in a rise in motor temperature and subsequent malfunction or fire.
- Some level of noise during usage is normal because of the tapping function.
- Stop operation immediately if you experience any discomfort.
- This product is not intended for self-diagnosis or self-treatment of conditions that should be managed by a qualified health care provider.
- People under medication or with medical conditions, please consult your doctor before use. Please do not use the product if you are not feeling well or if you are pregnant.


The producer is not responsible for damages caused by wrong use of device.

# INSIDE THE PACKAGE

The package includes: iNeck3 Massager, Power Adapter, User Manual, Fabric Sleeve.

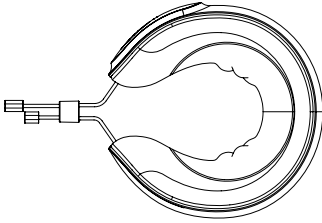iNeck3  
Massager

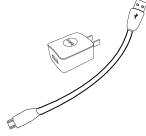Power Adapter


Fabric Sleeve

User Manual

# FUNCTION DEPLOYMENT
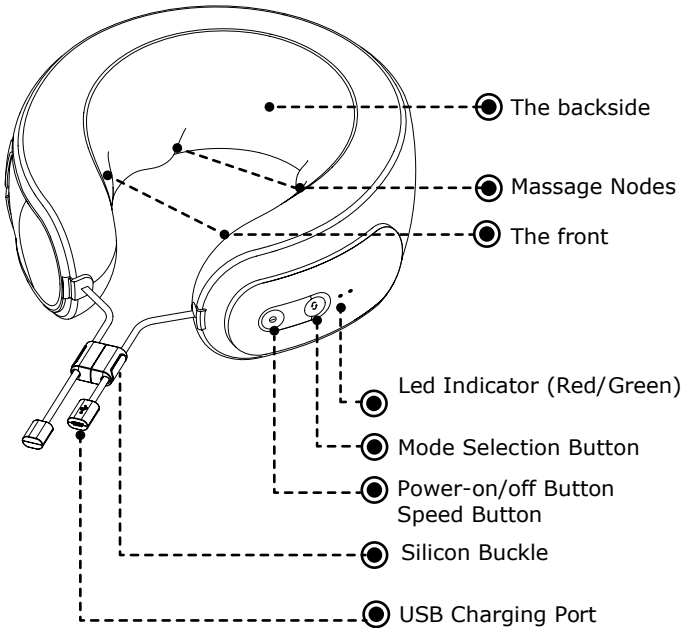
# PREPARATION

Prior to use, check the condition of the massager and accessories.

1. No damage on power charger and the massager.
2. Fully charge the battery prior to use.

## Charging

1. When the power is low, insert the USB cable to the USB charging port in the strap. Connect to the adapter to the power outlet.
2. The Battery Power Indicator flashes RED then stops flashing when fully charged.
3. Charging Time is about 10 hours.
4. Please do not use the massager during charging.

## Wear

1. Pull the rope out of the Silicon Buckle.
2. Place the main unit on the cervical vertebra (back of the neck).
3. Adjust massager to the most comfortable fit and insert the rope to the Silicon Buckle.

## Power-on

Press the  button. The LED light turns green and the massager immediately starts.

(Preset default: Shiatsu Mode1 for 10 minutes, slow speed)

One beep for slow speed, two beeps for medium speed and three beeps for high speed. To switch among massage modes, press the mode button. The modes are listed in the mode section of this manual.

## Power-off

Press and hold the button . The LED light and massager turn off.

# OPERATION

Mode Selection

Note: Default Mode is Shiatsu **Mode1**

Press the  button the massage mode changes in this sequence:  
Mode2 > Mode3 > Mode1

The three modes differ in direction of strokes. See below:

**Mode1** the massage nodes keep rotating back and forth.

**Mode2** the massage nodes do not rotate during transition from the backside to the front side and rotate during transition from the front side to the backside.

**Mode 3** the massage nodes rotate during transition from the backside to the front side and do not rotate during transition from the front side to the back side.

Note : Pressing the button  continuously to stop the translation of massage nodes and select a fixed-position massage. After hearing a " beep " sound , keep pressing the button and select the position you prefer , then you release the button to enjoy the massage.

### **Automatic Shut-Off**

The unit will automatically shut o after 10 minutes.

### **Cleaning and Care**

To clean, remove the fabric sleeve and wash.

## **PRECAUTIONS**

1. Please read the manual carefully and operate accordingly.
- 2.This unit is not a toy. Keep away from children.
3. Keep the massager away from humidity, heat and anything flammable.
4. If your hair is long, we recommend that you tie it up before using the neck massager. This will prevent the hair from possibly getting caught in the machine
5. Use only the supplied power charger, use it with caution, and according to specifications.
6. This massager should not be used by:
  - a. People who suer from swelling, inflammation or blood stasis in the neck, or have undergone neck surgery.
  - b. People who suer from slow movement, bad physical condition or intellectual disability.
  - c. People who has stroke, hypertension or heart disease.
  - d. People who suers from long illness, infirmity and age.
  - e. People who has open wound, serious infection, skin disease or broken skin, blood disease or bleeding tendency.
  - f. Gravida and child.
7. Do not use the massager when you are on an empty stomach or full stomach or after drunk or after vigorous exercise.
8. The massager is a personal care product and does not substitute health care professionals' advice and treatment. In case of doubt or if necessary, please consult your health care professional prior to use.
9. Stop using the massager immediately if feeling uncomfortable.
- 10.Do not us the massager when driving.
- 11.Do not use it wherever it is high-temperature or wet.
- 12.To clean,take o the sleeve to wash and wipe the main unit with alcohol but not water.

13. Remove all neckwear, such as necklace , scarve etc.
14. Do not use the massager for pregnant women, children, menstrual women.

## TROUBLESHOOTING & TECHNICAL SPECIFICATIONS

When operational problems occur, please refer to the chart below. Do not disassemble the machine or try to repair it yourself.

Problem	Possible Reason	Possible Solution
Unable to boot	Lithium battery power is low.	Charge the unit.
Massage Strength or Speed is low	Lithium battery power is low.	Charge the unit.

### Technical specifications

Power Supply		DC5V --- 1.0A
Ideal Environment	Temperature	5℃ - 40℃
	Relative Humidity	≤80%
	Atmospheric Pressure	860hpa-1060hpa
Dimensions		Main Unit: 255mm(L)*265mm(W)*90mm(H)
Net Weight		990g

## STORAGE & TRANSPORTING

### Storage

1. If left unused for long periods of time, turn the power on every 6 months to keep the battery from dying.
2. Maintain good ventilation and keep away sunlight and flame.

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### **Protect the environment!**

Do not throw electrical items in the bin.  
Take them to the nearest recycling centre so that they can  
be destroyed or recycled.

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