**NATURE &
DECOUVERTES**

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
OCTOPUS MINI-MASSAGER

MINI-MASSAGEGERÄT OKTOPUS

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MINI-MASSAJADOR OCTOPUS

Réf. 15210380 - 15210390 - 15210400

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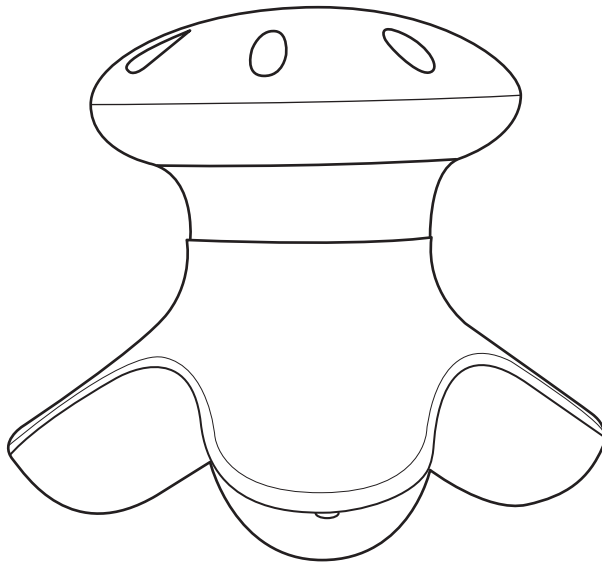
SCHEMA DETAILLE DES COMPOSANTS

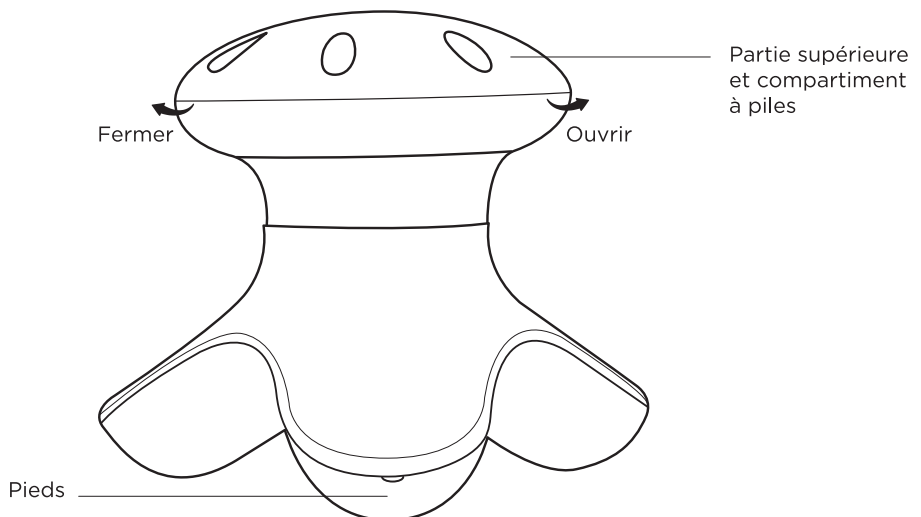PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

WARNING

- Do not place or store this appliance where it can be dropped into a tub or sink.
- Use this appliance only for its intended use as described in this manual.
- Never operate this appliance if it is damaged or it is not operating correctly.
- Close supervision is necessary when this appliance is used by children, invalids or disabled persons. This device should not be used by children without adult supervision.
- Turn OFF appliance when not in use.
- Do not place or drop into water or any other liquid.
- Keep appliance away from heated surfaces.
- Do not use this massager as a substitute for medical attention.
- This appliance is designed for household use only.
- If you have concerns regarding your health, consult a doctor before using this product.
- Do not mix old and new batteries.
- Individuals with pacemakers should consult a doctor before use.
- Those with the following conditions should consult their doctor: Pregnancy, Diabetes, Heart Diseases, Malignant Tumor, Conditions requiring rest.
- Never use this product directly on swollen or inflamed areas, skin eruptions or an open wound.
- If at any time an increase in pain or rash occurs, discontinue use and consult your physician.

The producer is not responsible for damages caused by wrong use of device.

WORKING

1. Unscrew the colored cap counterclockwise.
2. Insert the 3 AAA batteries into the compartment and put the cover back.
3. Simply turn the mini massager on with some pressure on the top
4. You can now use the mini massager on your whole body thanks to the soft vibrations: neck, lower back, shoulders, legs, feet and ankles. It can be used in the shower too to combine the benefits of massage and warm water.
5. After completing your massage, release the pressure on the top, it will automatically turn the unit off.

Remark: A massage session longer than 15 minutes is not recommended.


MAINTENANCE

- Keep it clean when you finish to use. You can clean the massager using a soft damp cloth to wipe the device to remove any dirt. If needed use some alcohol.
- Do not expose to moisture and strong sunshine.
- Take the battery out when you plan to stop using it for a long time.
- Do not put heavy weight on the machine.
- Do not disassemble the instrument.

TECHNICAL SPECIFICATIONS

- Model: Mini Massager
- Output power supply: 4.5V 0.29A
- Weight: 205g
- Batteries: 3 x AAA batteries (not included)
- Dimensions (approx) : H10.5 x 9.5 cm

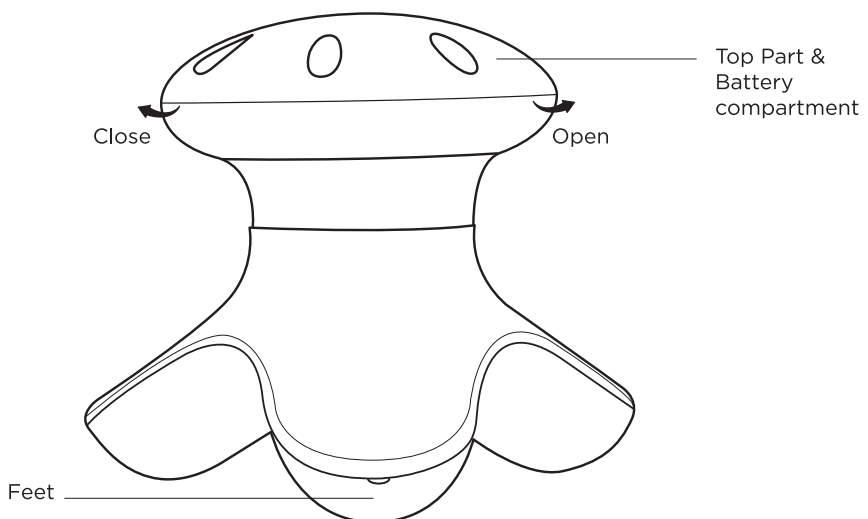DETAILED DIAGRAM OF THE COMPONENTS

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