NATURE &
DECOUVERTES

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PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

	Do not use folded or rucked
	Do not insert pins

- When the appliance is likely to be used for an extended period, for instance by the user falling asleep, controls should be adjusted to a setting recommended for continuous use.
- Examine the appliance frequently for signs of wear or damage. If there are such signs, if the appliance has been misused or does not work, return it to the supplier before switching it on again
- Do not use if wet

SAFETY GUIDELINES

- This appliance must not be used by persons insensitive to heat or by other vulnerable persons incapable of reacting in the event of overheating.
- Children under three years of age must not use this appliance due to their inability to react in the event of overheating.
- This appliance can be used by children aged eight or over and by persons with reduced physical, sensory or mental capacities or those lacking experience or knowledge, provided that they are properly supervised or the instructions for safe use of the product have been provided to them and the risks involved have been understood. Children must not play with the appliance.

	Not to be used by very young children (0-3 years)
	Read the instructions

- The appliance can be used by children older than 3 years and younger than 8 years under supervision and with the control always set to minimum temperature value.
- The cable and control unit of the appliance can give rise to risks of entanglement, strangulation, tripping or treading if not correctly arranged. The user shall make sure that excess ties and electric cords shall be arranged in a safe way.

- The appliance must not be cleaned or maintained by unsupervised children.
- The appliance should not be used by children or persons with disabilities or when asleep or by persons insensitive to heat.
- Do not fall asleep while using this appliance.
- Do not use on areas of the body which are swollen, injured or irritated.
- If you experience pain or discomfort while using the appliance, stop using it immediately. Not intended for medical use in a hospital.
- This appliance is not a toy, please ensure children do not play with it.


- This appliance operates via an electrical power source. Therefore, please observe the safety regulations for all electrical appliances.
- Only connect the appliance to the indicated voltage.
- Only use with the control unit supplied (model GL-01).
- To ensure the proper functioning of the control unit (thermostat), do not place it on top or underneath the appliance.
- Do not use it if its creased or crumpled.
- Before turning on the power, make sure that all components of the appliance are in perfect working order. In the event of a fault or malfunction, immediately turn off the power and unplug the appliance.
- When the appliance is not being used, store it away as instructed.
- Do not crease the appliance by placing objects on it during storage.
- Do not crumple the appliance. When storing, fold it carefully.
- When not in use, store the appliance in a dry place.
- When storing the appliance, allow it to cool down before folding it.
- If the power cable is damaged, it must be replaced by the supplier, one of its maintenance agents or any other qualified person in order to avoid any danger.
- Do not insert pins or pointy objects into the appliance.

USE

Thank you for choosing one of our products. We hope this product will bring you immense satisfaction.

- This appliance is for personal use only. It cannot be used for commercial purposes, in hospitals or medical practices.
- If you suffer from a pathology (varicose veins, thrombosis, etc.), please check with a doctor before using this appliance.
- Do not use this appliance if you have open wounds, burns, bruises, swelling or other similar problems.
- If you experience any discomfort, stop using the appliance. Heat-sensitive persons must use this appliance with caution and stop using it when they feel uncomfortable.
- Never use this appliance on animals.
- Only use this appliance in a dry room and on a dry body. Do not use in the bath, shower, above a sink full of water, in a swimming pool, sauna or outside.
- To benefit fully from the lumbar heat belt, please read this instruction manual carefully and follow the instructions regarding maintenance and safety. In this way, your appliance will last longer.
- Keep this instruction manual in case you need it for future reference, or if you give the appliance to another person.

CONTROL UNIT

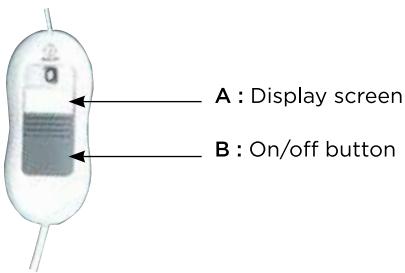A : Display screen
B : On/off button

1. Electrical connection:

To switch on the heat belt, plug the appliance plug into the mains socket (only plug into an appropriate outlet) and turn on the power. Please ensure that the voltage of the appliance (see data plate) and the mains voltage correspond.

2. Use:

Preferably apply the heat belt to the following area: the back.

The LED light will blink when the auto-off time of 90min is reached. The end user can reset the product by switching the controller to "0", and then turn it to level "1", "2" or "3".

CAUTION: Always use the heat belt with its cover. It provides a minimum layer to avoid overheating. Do not cover the remote control. This may lead to malfunctions.

REGULAR STORAGE AND INSPECTION

1. Unplug the remote control from the socket.

CAUTION!

- Store the appliance in a dry place.
- Do not place objects on the heat belt when storing it in order to avoid creating large creases.
- Regularly check the appliance and the power cord for signs of wear and tear or damage.
- In case of wear and tear or malfunction, please take the unit to the sales outlet.

CLEANING AND MAINTENANCE

The heat belt is machine washable at 30°C on a gentle cycle.

- Do not bleach.
- Do not put the heat belt in the dryer.
- Do not iron the heat belt.

3. Power supply:

Turn the switch on the remote control to the desired setting. The belt will turn on and begin to heat up.

NOTE: When using the appliance for the first few times, it may emit an unusual odour. It is harmless and will disappear during use.

4. Switch on the remote control:

You can select the follow settings:

Setting	Use
0	Stop
1	Low heat
2	Medium heat
3	High heat

NOTE: you can check the setting in the control window.

5. Turning off the appliance:

Turn the switch on the remote control to 0 (stop). Then remove the plug from the socket.

Made of fabric:

100% polyester, 100% polyester filling.

HEAT BELT:

1. First unplug the belt plug from the power supply.

2. You can machine wash the belt at 30 °C on a gentle cycle. Use a washing detergent for delicate laundry and use the dose recommended by the supplier.

• The belt should not be cleaned with chemicals, dried, machine dried or ironed.

• Please note, it is not recommended to wash the heat belt too frequently. Therefore, you should only wash the belt in the washing machine a maximum of 5 times during its lifespan. We recommend that you clean the unit by hand as much as possible.

3. Just after washing, stretch the wet belt to its original size and leave it dry out flat on a dryer.

CAUTION:

- Do not use clothes pegs or similar items to hang the heat belt on the dryer - they may cause damage to the belt.

CAUTION:

- Do not reconnect the remote control to the heat belt until the plug and the belt are completely dry - it cause damage to the belt.
- Never turn on the heat belt to dry it. You risk getting an electric shock.

REMOTE CONTROL:

Clean using a dry cloth only.


TECHNICAL SPECIFICATIONS

Model	HP-007A	Weight	0,44 kg
Nominal power	100 W	Nominal voltage	220-240V / 50Hz
	For indoor use only		This is a double-insulated appliance, which corresponds to Class 2 protection.
	When this appliance is at the end of its lifespan, please dispose of it in accordance with the Waste Electrical and Electronic Equipment Directive 2002/96/EC (WEEE).		

Subject to technical and design modifications as part of continuous product development.

NOTE REGARDING COMPLIANCE WITH DIRECTIVES

The manufacturer hereby declares that the HP-007A is in compliance with the main requirements of the European Electromagnetic Compatibility Directive (2014/30/EU) and the Low Voltage Directive (2014/35/EU).

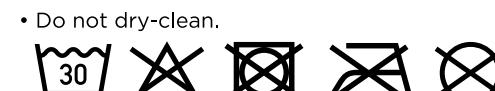Follow the cleaning instructions on the label.

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