**NATURE &
DECOUVERTES**

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PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

AVERTISSEMENT

The lightning flash with arrow-head symbol within an equilateral triangle is intended to alert the user to the presence of uninsulated “dangerous voltage” within the unit’s enclosure that may be of sufficient magnitude to constitute a risk of electric shock.

The exclamation point within an equilateral triangle is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the unit.

- Use this massaging mattress only for its intended use as described in this manual.
- Do not expose it to humidity, dust, or extreme temperatures.
- Always unplug this massaging mattress from the electrical outlet immediately after use and before cleaning.
- Do not place, drop or submerge in water or other liquid.
- Make sure that any objects do not fall and liquids are not spilled onto the unit.
- Never use pins or other metallic fasteners with this appliance.
- Never drop or insert an object into any opening.
- Carefully examine the cover before each use. Do not use the appliance if there are any signs of deterioration such as cracking or blistering.
- Keep the unit dry—do not operate in a wet or moist environment.
- Do not use wet hands to connect or disconnect the power plug from the electrical outlet.

- For adult use only.
- Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled individuals.
- An appliance should never be left unattended while plugged in. Unplug when not in use.
- Consult your physician before using this massaging mattress.
- Do not use the product if you are not feeling well.
- This product is not intended for self-treatment of conditions that should be managed by a qualified health care provider.
- The massaging mattress is designed as a non-professional massager to soothe tired and aching muscles. It is not intended to be a substitute for medical attention. If your ailments and symptoms persist or worsen after use, contact your physician immediately.
- Do not use the massaging mattress in combination with other therapeutic equipment or electric blanket, etc., as this may result in ineffectiveness or injury.
- Unplug the massaging mattress during lightning storms.
- Protect the power cord from being walked on or pinched.
- Keep cord away from heated surfaces. Wrap the cord loosely around the unit when storing.
- Do not operate under a blanket or pillow. Excessive heating can occur and cause electric shock or injury to persons.
- **IMPORTANT:** Any individual who may be in the following

- cases should consult a medical professional before using this massaging device designed for home use: be pregnant or recently gave birth, has a heart condition or a pacemaker, suffers from diabetes or other sensory impairment, phlebitis and/or thrombosis, joint dysfunction, inflammation or injuries, weak bones, osteoporosis, a high fever, an abnormal or curved spine, is at an increased risk of developing blood clots, has pins/screws/artificial joints or other medical devices implanted in his/her body, is being treated for a medical condition or physical injury or is under a doctor’s care.
- Do not use the massaging mattress if wearing loose clothing or jewelry.
- Keep long hair away from the massaging mattress while in use.
- Do not use on an infant or invalid or on a sleeping or unconscious person.
- Do not use on insensitive skin or on a person with poor blood circulation.
- The unit should be situated away from direct sunlight or heat sources such as radiators, electric heaters, stoves, or other units.
- Never block the air openings of the unit with materials such as clothing, plastic bags or papers, or place it on a soft surface such as a bed, where the air openings may be blocked. Keep the air openings free of hair and other things.
- Do not overload electrical outlet. Use only the power source as indicated.
- Do not carry this unit by cord and do not use the cord as handle.
- Never operate this unit if it has a damaged cord or plug, if it is not working properly, or if it has been dropped or damaged, or dropped into water. If the unit’s power supply cord or plug is damaged, do not attempt to fix it yourself.
- To avoid the risk of electric shock, do not disassemble or attempt to repair the unit. Incorrect repair can cause risk of electric shock or injury to persons when unit is used.
- This product is intended for household use only. Do not use this unit outdoors.

- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not operate in the presence of explosive and/or flammable fumes.
- Never remove the plug from the socket by pulling the power cord.
- To disconnect, ensure all controls are off and remove the plug from the electrical outlet.
- Do not stand on the massaging mattress.
- Do not insert your hand or foot along the paths of the massage rollers during use as this may cause injury.
- Stop operation immediately if you experience any discomfort.
- The product is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.
- Do not operate for more than 20 minutes continuously. It is recommended to turn off the massaging mattress and allow it to rest (cool) for 40 minutes before re-using. This will prolong the life of the unit.
- Do not perform continuous massage on the same spot of your body for more than 5 minutes at a time as this may result in excessive stimulation and may have adverse effects.
- Children should be supervised to ensure that they do not play with the product.
- Use the heating mode carefully. May cause serious burns. Do not use over sensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or disabled persons may be dangerous.
- Temperatures sufficiently high to cause burns may occur regardless of the control setting. Check the skin in contact with the heated area of the massager frequently to reduce the risk of blistering.
- Never use the massaging mattress on open skin eruptions, swollen or inflamed areas.
- Do not operate after taking pain killers, sedatives or alcohol.
- Do not apply massage rollers directly on the head, el-

- bow or knee joints, torso or abdomen.
- Burns may result from improper use.

- NEVER REMOVE THE COVER. Do not use the massaging mattress when the cover is removed. This may cause injury or malfunction.

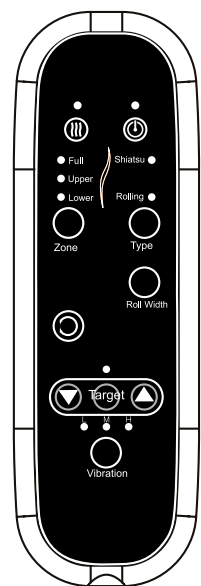
OPERATION

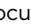
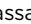
MASSAGING MATTRESS OPERATION

1. Place the mattress on a flat surface.
2. Insert the connector plug of the adapter into the socket of the massaging mattress.
3. Plug the AC adapter into an electric outlet.
4. Press the Power button to turn the massager on.
5. Press the Power button to turn the massager off. When the massager is powered off, all the LEDs on the remote will blink while the massager nodes move to the home position.

Note: the massaging mattress will automatically shut off after 20 minutes of operation. Do not operate for more than 20 minutes continuously

FUNCTION BUTTONS

1. **ZONE:** press the ZONE button to select: Full (entire upper body and back but not up to the lumbar), Upper (neck and shoulder blades) or Lower (waist and lower back but not up to the lumbar) massage coverage.
2. **MASSAGE NODE DIRECTION:** Press the  button to change the direction of the massage nodes when using Shiatsu massage.
3. **TARGET:** press the TARGET button to turn the target function on and off. Press and hold the  and  buttons to focus massage to a targeted area.
4. **HEAT:** For a soft feeling of warmth, press the  button to turn this option on and off. The heat can only be used when the massager is powered on. The temperature is 42 ° and cannot be adjusted. It works independently of the massage function. The user can therefore use the massage function without heating. It switches off automatically after 15 minutes.
5. **POWER:** press the  button to turn the massager on and off.
6. **VIBRATION:** press the VIBRATION button to select Low, Medium or High vibration massage. The vibration mode is only located on the second part of the mattress (the lower part of the body).
7. **ROLL WIDTH:** Press the ROLL WIDTH button to adjust the width of the roller massage.
8. **TYPE:** Press the TYPE button to select Shiatsu or Roller massage.

CARE AND MAINTENANCE

Composition: Polyurethane, Polyester and Foam.

- Washing not allowed.
- Bleaching not allowed.
- Do not dry in tumble dryers.
- Do not iron.
- No chemical cleaning.

Follow the cleaning instructions on the label.

CLEANING:

1. CAUTION: Unplug the mattress from the electrical outlet before cleaning.
2. Do not immerse the device in water.
3. Clean any stains with a damp cloth.

Cleaning and maintenance should not be done by unsupervised children.

STORAGE

1. Fold the mattress and store it with the rear strap.
2. Store in a cool, dry place.

SPECIFICATIONS

Power Supply	AC 120V-240V, DC 12V 4A
Power Output	42W
Power Frequency	50Hz
Weight	Approximately 8kg
Maximum user weight	110kg
Size	170 x 60 x 11 cm

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