
**NATURE &  
DECOUVERTES**

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**EYE MASSAGER  
MASAJEADOR DE OJOS  
OOGMASSAGEAPPARAAT  
MASSAJADOR DE OLHOS**

**Réf. 15218440**

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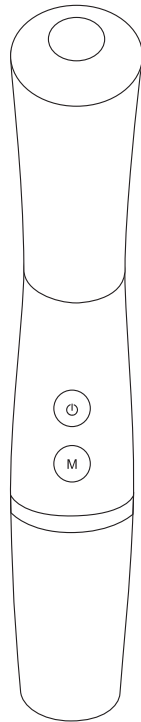# PLEASE READ THIS MANUAL CAREFULLY AND KEEP IT IN A SAFE PLACE


## WARNING

- Remove all packaging material from the appliance before use.
- This appliance shall not be used by children of 0 to 8-year old.
- This appliance can be used by children older than 8-year old and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, but only if they have been given supervision or instruction concerning use of the appliance in safe way and understand the hazards involved.
- Cleaning and maintenance must not be done by children.
- Have any repairs carried out solely by a qualified electrician. Never try to repair the appliance yourself.
- Do not handle the appliance with wet hands.
- Children shall not play with the appliance.
- Do not leave the appliance unattended when connected to the mains supply.
- Only use the attachments/accessories supplied with the product.
- Make sure that the appliance, the USB cable do not make contact with water.
- Make sure that the appliance, the USB cable do not make contact with hot surface, such as a hot hob or naked flame.
- Check the appliance's USB cable regularly to make sure it is not damaged.
- Do not use the appliance if falls into water.
- Never immerse the unit, the USB cable in water or any liquid.
- Do not place or store the unit near water, or any other liquid, or in a high humidity place.
- Do not use the liquids, such as water or flammable solvents to clean the product, please use a soft cloth to clean when product not in running.
- The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- The product must not be used by people who are sensitive to heat.
- Do not hand over this unit to the people or children who are not familiar using this product to avoid any risk of scalding.
- Continuous use of the appliance may, in the long term, damage the product and shorten its life span.
- Keep the appliance out of the reach of children and animals.

The manufacturer is not responsible for any damage caused by improper use.

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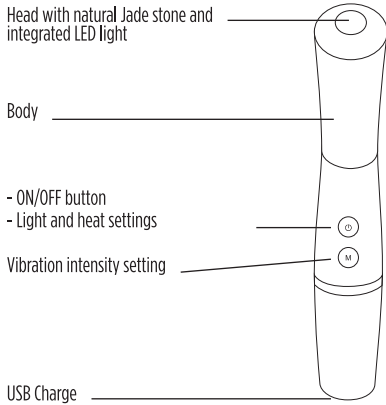- Operating time: 10 minutes before switching off automatically
- Voltage, current: 3.7V 30mA
- Power: 0,8 W
- Vibration frequency: 8000-12000 times/sec
- Temperature: 39°C-45°C
- Dimensions: 14,1 x 2,6 cm

## HEAT SETTINGS

BLUE LIGHT Low heat 39°C	RED LIGHT Medium heat 42°C	GREEN LIGHT High heat 45°C
-----------------------------	-------------------------------	-------------------------------

## OPERATION

1. Charge the Eye Massager  
Connect the USB cable to the unit and to a power source.  
When the battery is charging, the light on the Jade is RED.  
When the battery is ready for use, the light on the Jade is BLUE.
2. Apply your favorite cream or serum on your face.
3. Long press on the  button to turn on the device.
4. Short press on the  button to adjust the lights and temperature modes  
1<sup>st</sup> press: red light 42°C  
2<sup>nd</sup> press: green light 45°C  
3<sup>rd</sup> press: blue light 39°C
5. Press on  to change the vibration intensity from low to high
6. Massage the head over the cream or serum gently for 3-5 minutes until fully absorbed.
7. Long press on the  button to turn off the device.

## SAFETY MEASURES

The following safety measures must be observed when using this product:

1. Avoid any direct contact with the eye.
2. Prevent the light from shining directly in the eyes.
3. If you feel pain or a strange sensation on your skin, please stop using this unit immediately.
4. DO NOT let heat-sensitive people and children use this unit.
5. DO NOT wash with water.
6. This product is not recommended to be used by a pregnant woman or a person wearing a pacemaker.

## MAINTENANCE AND TROUBLESHOOTING

### Maintenance:

- Turn off the device after use.
- Wipe the head gently with a clean cloth after every use.

### Storage:

- Disconnect the adapter before storing.
- Store the device in a cool, dry place.

### Troubleshooting:

Problem	Possible cause	Solution
The device does not work	The power button is not on.	Turn on the power button with a long press on the power button.
	The device is out of battery	Charge the device
The device turned off	The auto shut off function of the device turns the device off after 10 mins.	Turn the device back on and be careful not to exceed 3-5 minutes of usage depending on your skin type.
The temperature is too high	The selected setting is too high for your type of skin	Choose a lower setting with a short press on the power button

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