NATURE &
DECOUVERTES

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REFLEXOLOGY MASSAGER

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PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

WARNING

- This appliance can only be used with the power cord provided.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- The power cord cannot be replaced. If the cord is damaged, please contact the seller to buy a new one.
- Failure to follow instructions may result in physical injury or material damage (electric shock, hurt, fire, etc.). These safety guidelines and warnings are intended not only to protect your health and that of others, but also to ensure the product is used appropriately. For this reason, please follow these safety guidelines and make them available to other persons when they use the product.
- This product must not be used by children, persons insensitive to heat or by other vulnerable persons incapable of reacting, due to their inability to react, in the event of overheating.
- If you experience pain or discomfort while using the product, stop using it immediately. Not intended for medical use in a hospital.
- This product is NOT A TOY. Children must not play with the appliance.
- The appliance must not be cleaned or maintained by unsupervised children.
- Use this appliance only for its intended use as described in this manual.
- This appliance is not intended for use as a medical purpose.
- Do not use this product if it is wet.
- The appliance should never be left unattended when plugged in. When not in use, unplug the appliance from an electrical outlet and store the appliance in a dry place.
- After use, please unplug the power plug.
- Connect the appliance only to the socket with appropriate voltage.
- The use of power supply with the voltage outside the specified range may cause accidents such as electric shock, fire, etc. Always use power supply with voltage within the specified range.
- Before turning on the power, make sure the machine is on a flat surface and not tilted.
- Never use this appliance and unplug immediately if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.

- Do not carry this appliance by supply cord or use cord as a handle.
- Keep out of the reach of children.
- Keep sharp objects away from the appliance to avoid damaging it.
- Protect the device from heat, water, naked flames, direct sunlight, extreme temperatures, humidity and collisions.
- Don't put hands into the machine to use, it is intended for foot use only.

It is forbidden to put hands into the cloth sheath

-  If the pins of the plug parts are damaged, the plug-in power supply shall be scrapped.


These instructions are also available on Nature et Découvertes website: www.natureetdecouvertes.com


SAFETY PRECAUTIONS

- Please use the machine for not more than 30 minutes, otherwise it may cause pain to muscle and nerve.
- The following users should consult a physician before using the product:
 - A user who has a pacemaker or similar implantable medical device;
 - A user who suffers from heart disease;
 - A user who suffers from muscle strain;
 - A user with fragile or injured bones, especially a user with injured back bone;
 - A user who just underwent surgery;
 - A user who is pregnant or just gave birth;
 - A user who suffers from osteoporosis;
 - A user who suffers from blood circulation problems or thrombosis;
 - A user who suffers from skin infection;
 - A user who receives treatment or adjuvant treatment;
 - A user with abnormalities in the limbs;
 - A user who suffers from spondylosis or vertebral injury caused by accident or disease.
- Do not use this product in extremely humid and high-temperature rooms nor use it beside heating objects (such as electric heating, stoves, etc.).
- Before or during the use of this product, please do not apply medical or cosmetic creams or lotions to your feet.
- Do not stand on the product or put a heavy object on it.
- This product is not for use in a vehicle.
- Electrical appliances are dangerous when they get in contact with water. Please do not use this product in a humid environment like bathroom or places near a faucet or swimming pool.
- If the power cord is damaged, it must be replaced by the supplier or its service center or a person with relevant qualifications, so as to prevent risks.
- This product contains a lithium button cell battery. A lithium button cell battery can cause severe internal chemical burns if swallowed. Discard used batteries immediately. Keep new and used batteries out of the reach of children. If batteries have been swallowed or placed inside any part of the body, seek medical advice immediately.

- To change the battery: using a screwdriver, unscrew the screw to open the cover. Remove the old button battery. Insert a new button cell battery into the battery compartment, ensuring correct polarity.
- Button batteries, like any battery, must be disposed of at a dedicated collection point for recycling. The batteries must not be recharged. Do not combine different types of batteries or new and used batteries. Pay attention to the polarity (+ and -) when installing the batteries. Used batteries must be removed from the product. The power terminals (the battery) must not be short-circuited.
- The remote control is powered by a CR-2025 button cell battery (battery not included).

PRODUCT VIEW

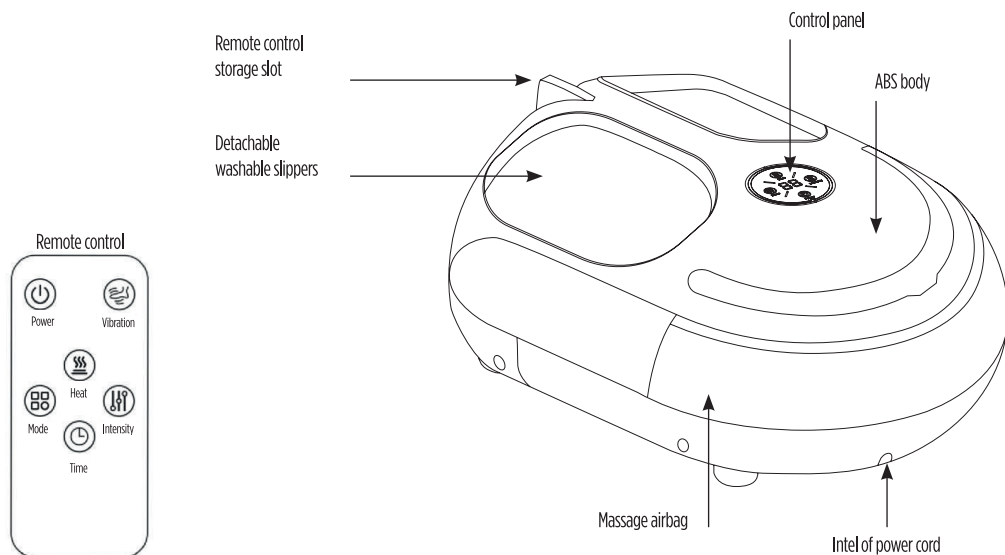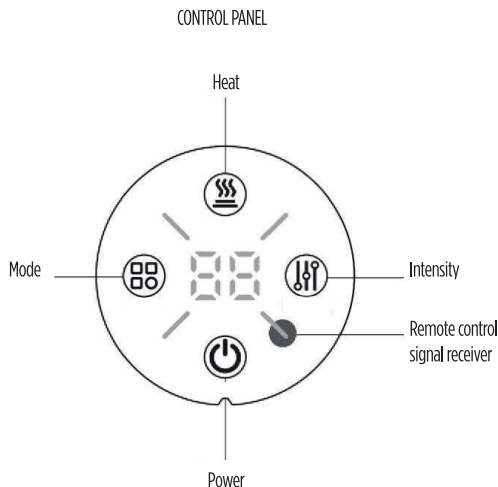Long press 2 seconds to start the machine, then short press to adjust the timing (5, 10, 15, 20, 25 or 30 minutes). Default timing is 15 minutes.

Mode button: there are 4 different modes to choose from. Short press to select between mode 1, 2, 3 and 4.

Mode 1: forward rolling massage + air pressure.

Mode 2: reverse rolling massage + air pressure.

Mode 3: forward and reverse rolling massage only.

Mode 4: air pressure only. Default mode is 1.

(Optional): long press to select vibration mode (on/off). Vibration is off by default.

Intensity button: 5 levels of massage pressure. Short press to select the desired pressure. The lowest pressure is set by default.

2 levels of gentle heat: short press to select the desired heat. Heat is off by default.

This appliance has a timer. When the timer expires, the appliance automatically switches to standby mode. Its energy consumption then drops to 0.4 W.

TECHNICAL DATA

- Composition:
removable slippers: 100% polyester,
massager: ABS
- Voltage: 100-240 V, 50/60 Hz
- Power: 40 W
- Timer: by default, the massager is set to 15 minutes.
You can change this to 20, 25 or 30 minutes, as you wish.

The company reserves the right to improve the product and the manual. If there is any change without notice, the appearance of the product in kind shall prevail. This book has been carefully checked, but there may be printing mistakes or misunderstanding of the content, please understand. The company reserves the right to explain.

MAINTENANCE

- If there are stains on this machine after use, please wipe it with a cloth soaked in water diluted neutral detergent. Do not use chemicals such as petroleum essence, alcohol or diluent on the massager. Then, let it air dry. Please do not dry the machine in an intense manner.
- Warning: before maintenance and cleaning, please turn off the power switch and unplug the power plug.

CLEANING

- Care instructions: disconnect from the power supply before cleaning. Unzip and remove the removable slippers to wash them. Hand wash only.

TROUBLESHOOTING

Troubles	Fix
The machine is not working after the connection to the power supply.	Check that the appliance power cord is not damaged.
	Check if the power cord is correctly plugged into the machine, and if the plug of the power cord is correctly plugged into the power supply.
	Check that the product and motor temperature are not excessively high after a long period of operation. If this is the case, the temperature safety device will be triggered automatically. Switch off the appliance and allow the motor to dissipate the excessive heat for at least 30 to 50 minutes.
The control panel is unresponsive and the machine has stopped working.	Please press the power button to turn off the machine and then restart again.
The remote control does not work.	Replace the CR-2025 button battery (battery not included).

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