
**NATURE &
DECOUVERTES**

MASSEUR PÉTRISSAGE TRAPÈZES
TRAPEZIUS PETRISSAGE MASSAGE DEVICE
MASAJEADOR DE TRAPÉCIOS CON
TÉCNICA DE AMASAMIENTO
TRAPEZIUM KNEED-MASSAGE-APPARAAT
APARELHO DE MASSAGEM PARA
MASSAJAR OS TRAPÉZIOS

Réf. 15235590

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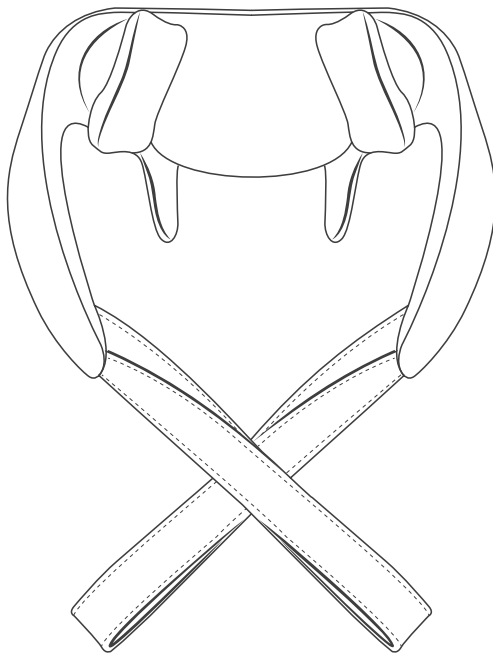PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

This set contains:

- 1 trapezius petrissage massage device
- 1 mains adaptor:
input: AC 100-240 V 50/60 Hz,
output: 5.0 V $\equiv$ 2.0 A 10 W
- Instruction for use
These instructions are also available on
Nature et Découvertes website:
www.natureetdecouvertes.com

INSTRUCTIONS BEFORE USE

- This appliance can be used by children aged from 8 years and above and persons with reduced physical sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Children under the age of three are not to use this appliance due to their inability to react to overheating.
- Do not crease the appliance by placing items on top of it during storage.
- If the power cable is damaged, it must be replaced by the supplier, an approved repair center or any other qualified person in order to avoid any danger.
- Do not use the kneading massager if the fabric cover is damaged.
- Do not insert any object or flammable material into the kneading massager, this may lead to an electric shock, breakdown or risk of fire.
- Do not use excessive pressure during use of the kneading massager.
- Do not cover the kneading massager during use. This may cause the motor to increase in temperature and result in breakdown or fire.
- It is normal for a certain level of noise to be heard during operation. Immediately stop the appliance if you experience any discomfort.
- This appliance is not intended for self-diagnosis or the self-treatment of diseases that should be treated by a qualified therapist.
- Anyone taking medication or suffering from a medical disorder must consult their doctor prior to using the appliance. Do not use this product if you feel unwell or if you are pregnant.
- Unsuitable for use by young children 0-3 years.
- Read the instructions.
- When the appliance is likely to be used for an extended period, the settings should be adjusted to the recommended level of continuous use.
- Do not use wet.

- The external flexible cable or cord of this adaptor cannot be replaced; if the cord is damaged, the adaptor shall be scrapped.
- This appliance contains batteries that are only replaceable by skilled persons.
-  If the pins of the plug parts are damaged, the plug-in power supply shall be scrapped.
- The appliance is only to be used with the adaptor / power supply unit provided with the appliance.
- Do not use this product in close proximity to loose clothing or jewelry, keep long hair away from massager while in use.
- When not in use, store the appliance in a dry place.
- The appliance cannot work when the charging plug is inserted.
- Keep the appliance away from sharp objects or sharp edges at all times.
- Examine the appliance frequently for signs of wear or damage. If there are such signs, if the appliance has been misused or does not work, return it to the supplier before switching it on again.
- This appliance is not intended for medical use in hospitals.
- This appliance must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
- **WARNING:** for the purposes of recharging the battery, only use the detachable supply unit XSC-0502000SEUU provided with this appliance

Not intended for children ages 3 and under

Do not fold

Read instruction manual

No needling

FUNCTION BUTTONS

On/Off

Press the on/off button for 2 seconds to turn on and off.

Speed

2 speeds to choose from when the appliance is turned on, low speed is on by default, single short press to turn to high speed.

Mode





To adjust the kneading mode, 2 modes to choose from: continuous kneading and intermittent kneading.

Soothing heat Mode

Single short press to turn on and off the soothing heat function.

SUGGESTED USE

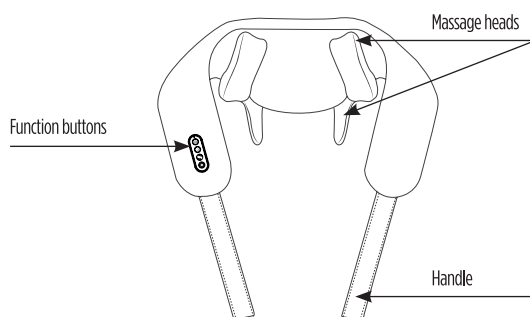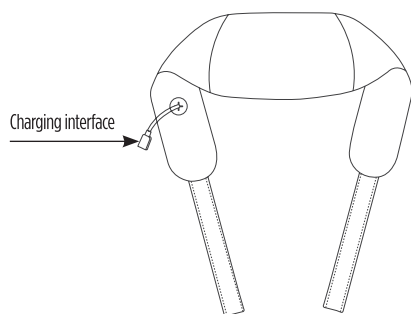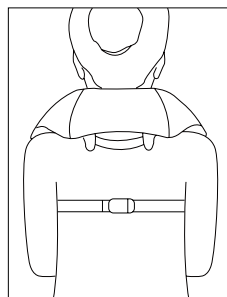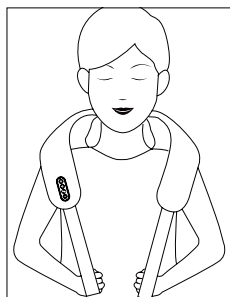Warning: when the product is in use, it needs to be pulled or held by hand only.

REMINDER

- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- The appliance is only to be used with the adaptor/power supply unit provided with the appliance.
- It must only be supplied at safe extra low voltage corresponding to the marking on the appliance.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

TROUBLESHOOTING

Trouble	Possible cause	Solution
The device does not work.	Batteries are out of power.	Charging.
	The ON/OFF button was not pressed correctly.	Press the ON/OFF button again.
	The device has switched off automatically after operating for 10 minutes.	Press the ON/OFF button again.

CARE AND MAINTENANCE

- Composition: synthetic leather 100% polyester coated with polyurethane.
- Do not wash.
- Do not iron.
- Do not bleach.
- Do not dry clean.
- Do not tumble dry.

- Massage system's protective fabric is removable and washable.
- Fabric: 77% nylon, 23% elastane.

TECHNICAL SPECIFICATIONS

- Item: trapezius petrissage massage device
- Adaptor model: # XSC-0502000SEUU
- Input Voltage(V): 100-240 V
- Input AC Frequency (Hz): 50/60 Hz
- Output voltage(V) DC: 5 V
- Output Current:(A): 2 A
- Power(W):10 W
- Autonomy : 50 minutes
- Works with rechargeable battery 7.4V 2000 mAh
- Average Active Efficiency (%): 79.89 %
- Efficiency at Low Load (10%)(%): ≥ 75.19 %
- No-load Power Consumption(W): ≤ 0.1 W
- Automatic timer: approximately 10 minutes
- Gross weight: 1.3 kg
- Net weight: 1.22 kg

BATTERY WARNINGS

- The product charging slot should not be short-circuited.
- Attention should be drawn to the environmental aspects of battery disposal, Don't throw used batteries in dust bin. Please dispose of it in accordance with the requirements of the local recycling authority.
- The batteries (built-in batteries) shall not be exposed to excessive heat such as sunshine, fire or the like.
- If battery leakage occurs, avoid contact with skin, isolate the leaking battery in a sealed plastic bag and put it in the scrap following the instructions regarding environmental protection.
- If battery leakage contact with the skin, mucous membranes or eyes, rinse thoroughly with water and contact your doctor or ophthalmologist.
- The equipment does not have a standby mode or other equivalent state in terms of energy efficiency requirements, power management function.
- The power consumption of the equipment is likely to be higher than other equipment models meeting these functional requirements.

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