NATURE &  
DECOUVERTES

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**EN**

**PLEASE CAREFULLY READ  
THIS MANUAL AND KEEP IT  
IN A SAFE PLACE**

**PRECAUTIONS FOR USE**

- This appliance shall not be used by children.
- Cleaning and user maintenance shall not be made by children. Appliances can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and if they understand the hazards involved.
- Children shall not play with the appliance.
- Always disconnect the appliance from the supply before assembling, disassembling or cleaning.

- Inspect the appliance frequently for signs of wear or damage. If any such signs are found, if the appliance has been used abnormally or if it does not work, return it to the supplier before using it again.
- Do not use this appliance near bathtubs, showers or water-filled reservoirs.
- Use of this appliance may cause discomfort or adverse effect to, persons suffering from skin disease or inflammation of the skin, persons with sensitive skin and persons having an allergic constitution. Stop using it if any discomfort occurs;
- This appliance is not to be used on the parts where the integrity of the skin is damaged (such as burns, blisters, lesions, incision, wounds, active skin disease or infection);
- This appliance shall not be used by, or on, persons wearing implants not

compatible with its use.

- This appliance shall not be used by or on persons up to and including 14 years of age.

WARNING: adverse effects may occur if the instructions are not followed;

- Prolonged or repeated application on one site may potentially damage the skin;
- Do not bypass or remove the inherent safety device of the appliance.
- Do not use the appliance in conjunction with flammable based cosmetic products, unless they are recommended by the manufacturer;
- Do not use the appliance on persons with an implanted electronic device (for example a cardiac pacemaker).
- Do not use the appliance close to the region of the heart.
- Do not treat certain parts of the human-

body such as eyes and genital areas.

- This device is not intended for medical use in hospitals.
- The appliance is only to be used with the power supply unit provided with the appliance.
- This appliance contains batteries that are nonreplaceable.
- For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.
- The appliance is only to be used with the power supply unit provided with the appliance.

GQ12-050200-ZG  
detachable supply unit

Read operator's manual

**WARNING:** Do not use this appliance near bathtubs, showers or water-filled reservoirs.

## WARNING

The following individuals should not use the device:

Failure to comply may cause accidents or other problems.

- If you have existing or suspected heart disease, or related disorders.
- Infants or if you are unable to express your intentions on your own.
- If you are unable to sense temperature or stimulation.
- If you have an autoimmune disease.

Consult a physician prior to use under the following conditions:

- If you have an acute disease.
- If you have an infectious disease.
- If you have a malignant tumor.
- If you have a febrile disease.
- If you have high blood pressure or a blood disease.

- If you have a fever.
- If you have sensory impairment due to severe peripheral circulatory dysfunction (e.g. diabetes).
- If you have a nerve damage.
- If you have a skin disease or atopic dermatitis.
- If you are a hospital outpatient and are taking medication.
- If you have a pre-existing skin condition, such as an allergic constitution or contact dermatitis, or if you have a concern about your skin when using this device.

Do not use in the following cases or on the following areas.

- If your skin is damaged, experiencing a rash, inflammation or other abnormalities.

- If you are under the influence of alcohol
- Head, anywhere inside the mouth, around the pubic region or any mucous membrane.
- Bulging moles.
- Viral warts.
- Areas where you have broken a bone.
- Areas of internal bleeding, swelling or scabbing after cosmetic medical procedures.
- Avoid using the device for 2 weeks after cosmetic medical procedures. After this period, consult with a physician and only use the device if the physician confirms there is no problem. (e.g. Hyaluronic acid injections, silicone injections, Botox, thread lifts, double eyelid surgery, HIFU, etc.).

- Areas where metal was implanted by surgery (reconstructive surgery, plastic surgery, cosmetic surgery).
- The power consumption was 0.1 W in the off mode.
- After the device is fully charged, the time after reaching the automatic shutdown mode is less than 1 min.


## PRODUCT INTRODUCTION

7 color Facial Wand is a beauty product that makes skin care by various technological such as LED light, high-frequency vibration, hot compress and EMS micro current.

## PRODUCT PARAMETERS

- Product name: 7-Color Facial Wand
- Product model: CCA-FM65A, 15235660
- Battery capacity: 3.7V 420mA
- Input voltage: 5 V
- Charging time: about 2 hours
- Product size: about 149 mm x 29 mm x 25 mm
- Net weight : about 45 g
- Power: about 0,9 W
- Maximum output: 18 V, 23,8uA

## GRAPHICAL DESCRIPTION
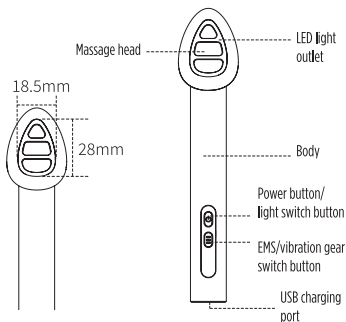
## OPERATING INSTRUCTION

1. Long press the power button to turn on the machine. The default mode is vibration + hot compress + LED red light.
2. When turning on the machine, the EMS and vibration modes are turned on at the same time. Short press the EMS/vibration button, the EMS microcurrent and vibration intensity will be switched from level 1 to level 2. Press it again to switch to level 3; 1-3 levels will be cycled.

3. Short press the power button to switch the LED light color. The switching sequence is red, green, blue, yellow, cyan, purple, white, and then cycle switch through the seven colors.
4. This unit automatically turns off after 10 minutes of care time.

## UTILISATION

1. Apply proper amount of cream or essence on the skin to keep the skin moist.
2. Select the desired LED mode and gear. For first-time users, it is recommended to use the vibration + EMS mode from the level 1.

### 3. Colours benefits

RED LIGHT: recommended for sensitive skin, helps soothe the skin

GREEN LIGHT : helps reduce skin redness and even out skin tone

BLUE LIGHT: recommended for sensitive skin, helps soothe the skin and reduce the appearance of pores

YELLOW LIGHT: helps reduce skin redness

PURPLE LIGHT : helps brighten the skin

CYAN LIGHT: helps improve the appearance of rough skin


WHITE LIGHT: recommended for skin with fine lines or wrinkles.

4. Gently press the massage head to improve the absorption of cream or essence, massage and lift the cheeks, forehead, chin, and other areas to improve your skin condition. stay at each point for 10 seconds, and repeat 2-3 times.

5. Other than facial skin, other body parts such as the neck, hands and legs can also be massaged.

Tips: the massage head can be rotated 180°, with the pointed head upward, it is suitable for the eyes; and the round head upward, suitable for larger areas such as cheeks.

6. After use, wipe the massage head with a piece of clean cloth, damp towel, or alcohol sponge.
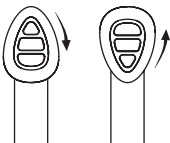
## ABOUT CHARGING

1. Plug the USB charging cable into the supplied 5V power adapter for charging.
2. When charging, the yellow light is on. When finished charging, yellow light goes off.
3. When in use, if the indicator light flashes frequently, it should be charged immediately as a power shortage reminder, otherwise it will stop working.
4. The charging cable should be properly kept and not lost.

## WARNING

1. This product is not waterproof. Please do not place this device in or near areas such as bathtubs, sinks, and wash-basins to prevent it from falling and damaging. To keep the product clean, please be sure to clean it after each use.
2. Ensure that the battery is fully charged at least once within 3 months, otherwise there is a risk of drying up

- and damage to the battery.
3. When we are using it, the massage head can be rotated, and the maximum rotation range is 180°.
  4. Do not use the device when charging.

## ATTENTIONS

1. If you have skin diseases, please consult a dermatologist before use.
2. For the first use, please try it on the arm.
3. If there is any discomfort, please stop using it immediately and consult a dermatologist.
4. For hygiene reasons, please do not share this product with others.
5. Place the device out of reach of children to prevent accidents.
6. Do not use corrosive chemicals to wipe the product.
7. Non professionals are not allowed to disassemble this product.
8. This product has a constant heating hot compress function, please use it with caution for temperature sensitive individuals.

## HOW TO CLEAN

To keep the product clean, please be sure to clean it after each use.

1. Disconnect the AC adapter from the main unit and the wall outlet.
2. To clean the device and electrode head, wipe them with a damp and well-wrung cloth to remove any dirt. Do not leave the device with lotion or gel on it.
3. Wipe the entire device with a dry cloth and store it after removing all moisture.

## COMMON PROBLEMS AND SOLUTIONS

This chapter provides an overview of the problems that may arise when using this product and proposes solutions. If you are unable to use the following information to solve any problems you encounter, please contact the customer service center in your country/region.

Question	Possible cause	Solution
The product cannot be used or suddenly stops working	Low battery	Try again after fully charging
	Product damage	Please consult after-sales customer service
The product overheats during use	Perhaps the ambient temperature is too high	Suggest standby and continue to use the device after the temperature returns to normal
After use, the massaged area turns red	Mild redness is a normal phenomenon that is harmless to the body and will quickly subside	Suspend use and wait for symptoms to subside before using the product again
Not very satisfied with the effect	Unable to feel the action of microcurrent on the skin	Skin is too dry, it is recommended to use with moisturizing products

When abnormal situations occur or faults cannot be eliminated, please contact the after-sales customer service in your corresponding language in a timely manner. Please keep the purchase invoice and warranty card for later use during product maintenance period.

## EXTERNAL POWER SUPPLY

- Manufacturer's name or trade mark, commercial registration number and address:  
Dong Guan City Gang, Qi Electronic Co., Ltd.,  
District, ShiYong, 532600, Heng Li Town, Dong Guan,  
Guangdong, PEOPLE'S, REPUBLIC OF CHINA
- Model identifier: GQ12-050200-ZG
- Input voltage: 100-240 V
- Input AC frequency: 50/60 Hz
- Output voltage: DC5,0 V
- Output current: 2 A
- Output power: 10.0 W
- Average active efficiency: 80.61%
- No-load power consumption (5V): 0.07 W

### Class II equipment

To identify equipment meeting the safety requirements specified for Class II equipment according to IEC 61140.

This symbol indicates that a safety connection to electrical earth (ground) is not required.

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