**NATURE &
DECOUVERTES**

LE COUSSIN CHAUFFANT POUR ÉPAULES
THE HEATING PAD FOR SHOULDERS
LA ALMOHADILLA TÉRMICA PARA HOMBROS
HET VERWARMINGSKUSSEN VOOR SCHOULDERS
A ALMOFADA TÉRMICA PARA OS OMBROS

Réf. 15236450

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IMPORTANT INSTRUCTIONS. RETAIN FOR FUTURE USE.

	Do not use folded or rucked		Not to be used by very young children (0-3 years)
	Do not insert pins		Read the instructions

- Examine the appliance frequently for signs of wear or damage. If there are such signs, if the appliance has been misused or does not work, return it to the supplier before switching it on again
- The appliance can be used by children older than 3 years and younger than 8 years under supervision and with the control always set to minimum temperature value.
- Recommended to select “Setting 1” for continuous use.
- The cable and control unit of the appliance can give rise to risks of entanglement, strangulation, tripping or treading if not correctly arranged. The user shall make sure that excess ties and electric cords shall be arranged in a safe way.
- When the appliance is likely to be used for an extended period, for instance by the user falling asleep, controls should be adjusted to a setting recommended for continuous use.
- Do not use if wet.

SAFETY INFORMATION

- Failure to follow instructions may result in physical injury or material damage (electric shock, skin burns, fire). These safety guidelines and warnings are intended not only to protect your health and that of others, but also to ensure the product is used appropriately. For this reason, please follow these safety guidelines and make them available to other persons when they use the product.
- This appliance must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
- Children under the age of three are not to use this appliance due to their inability to react to overheating.
- Do not use on areas of the body which are swollen, injured or irritated.
- If you experience pain or discomfort while using the product, stop using it immediately. This appliance is not intended for medical use in hospitals.
- This appliance is not intended for use in a hospital environment.
- Avoid direct and prolonged contact with skin. Do not place on bare skin.

- Do not insert pins or pointy objects into the appliance.
- Keep the appliance away from sharp objects or sharp edges at all times.
- Do not use the heating pad if it is pleated or creased.
- Do not use the heating pad if it is wet.
- When not in use, store the appliance in a dry place.
- Do not crease the product by placing objects on it during storage.
- When the appliance is likely to be used for an extended period the settings should be adjusted to the recommended level for continuous use.
- The appliance is only to be used with the detachable controller, that is marked on the appliance.
- Carefully check the heating pad before each use, if you notice wear or damage. Do not use if you notice wear, damage or signs of improper use on the heating pad, the switch or the cables. Return it to the supplier.
- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Do not fold or crease the heating pad during use. Neither pierce the pad with, nor attach safety pins or other spiky or sharp objects to the pad.
- Do not use the heating pad on children or persons who are disabled, asleep or insensitive to heat.
- Children should be supervised to ensure that they do not play with the appliance.
- Never use the pad unsupervised.
- Only use with the control unit supplied (model GL-20).
- The unit is designed for private home use only.
- Prolonged use of the heating pad at a high setting could lead to skin burns.
- Do not cover the pad with any other cushions. The control switch should not be covered or placed on or below the pad while the unit is in operation.
- Do not sit on the heating pad, only place the pad on the part of your neck you want to treat.
- Never touch a heating pad that has fallen into water. Unplug the unit from the mains outlet immediately.
- Keep the mains lead away from hot surfaces.
- Never carry, pull or turn the heating pad by the mains lead and do not let the lead get tangled.
- Do not use the heating pad when wet and only use it in dry environments (not in the bathroom or similar environments).
- The control and leads must not be exposed to any kind of moisture.
- If a fault occurs, do not attempt to repair the unit yourself. Repairs must only be carried out by an authorized specialist dealer or other suitably qualified personnel.
- If the power cable is damaged, it must be replaced by the supplier, one of its maintenance agents or any other qualified person in order to avoid any danger.
- When storing the heating pad, allow it to cool down before folding.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- For cleaning information and laundry symbols explanation, please see "CLEANING AND MAINTENANCE"

ITEMS SUPPLIED AND PACKAGING

The following parts are included as standard:

- 1 heating pad
- 1 instructions for use

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required.


If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

Please ensure that polybag packing is kept away from reach of children! Risk of suffocation!

OPERATION

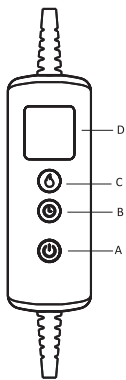A : ON/OFF

B : timer button

C : temperature button

D : digital display screen

1. Connect the AC power, the “” will be displayed on the screen which will disappear after 1-2 seconds.
2. Then the digital display screen will show “.” (Orange dot).
3. Press on/off button (button A), the digital display screen will appear “-”, the orange line will flash ten seconds. Press any button B or C except the on/off button within 10 seconds to switch on the product.
4. The heating blanket will be programmed to work under the default setting 3 (30 minutes) for timer and the 2nd heating level (temperature level 2) when using for the first time. The controller will remember the latest heating level and timer if the power plug has not been unplugged since the last time.
5. Press the button C to select your desired temperature from Level 1, 2, 3, 4, 5, 6 level, press button C to adjust the temperature level up and down. Press button B to adjust the timer level from 3 (30 minutes), 6 (60 minutes), 9 (90 minutes), H (120 minutes). Press B to adjust the timer up and timer down.
6. The LED display will indicate “O” when auto off.
7. Heating level 1 is recommended for continuous use.
8. Press button A to turn off the product when power on.
9. Repress the button A and repeat the step (3) to re-start the appliance when power off.
10. Please pull off the power plug when you are not using the product.

CLEANING AND MAINTENANCE

HEATING PAD:

- Before cleaning the pad, unplug it from the mains outlet and let it cool down for at least 10 minutes.
- The heating pad can be cleaned with a soft brush when dry.
- Never use strong detergents or hard brush.
- **Made of 100% polyester**
 - Washing, mild process, maximum 40° Celsius.
 - Do not bleach.
 - Do not tumble dry.
 - Do not iron the appliance.
- You should untangle the cable from time to time.
- When storing the appliance, allow it to cool down. Store the unit in a clean and dry place, laid out flat with nothing covering or resting on top of it.
- Do not dry clean.

Follow the cleaning instructions on the label.

SPECIFICATIONS

Name und Model: Heating pad 15236450

Power Supply: 220-240 V ~ / 50 Hz

Heat output: 120 W

Automatic switch off: 30/60/90/120 min

Switching level: 1-2-3-4-5-6

Protection class: 

Operating conditions: Only use in dry rooms as described in the operating instructions

Storage conditions: Laid out flat, clean and dry

Power consumption in off mode: 0,2 W.

The period after which the equipment reaches automatically off mode: less than 1 minute.

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Fabricant : DONGGUAN CITY GAOLAND TECHNOLOGY CO., LTD.