## **LE GUA SHA LED**

THE LED GUA SHA

EL GUA SHA LED

DE LED GUA SHA

O GUA SHA LED

**Réf. 15239210**

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# PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

## WARNING

- It is generally not advisable for people who have undergone facial plastic surgery, subcutaneous metal implant and those who are pregnant to use this product.
- Wipe the massage head with tissues or cloth after use. Do not wash it with water.
- Do not place the device in the bathroom or other humid environment, as this may cause a short circuit.
- DO NOT use alcohol, benzene or other solvents to clean the unit, as these substances may weaken and/or permanently damage the components and materials of the unit. If necessary, use only mild dishwashing soap and a damp cloth.
- Usage is prohibited under the following conditions:
  - Patients with malignant tumor
  - Patients with a bleeding tendency
  - Pregnant or postpartum women
  - Patients with acute suppurative inflammation
  - Users of heart pacemakers
  - Patients with metal implant
  - When applying the device near the heart
- This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- This appliance shall not be used by children. Cleaning and user maintenance shall not be made by children. Appliances can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and if they understand the hazards involved. Children shall not play with the appliance.
- Always disconnect and turn off the appliance from the supply before assembling, disassembling or cleaning.
- This appliance contains batteries that are non-replaceable.

Read instruction manual

WARNING: do not use this appliance near bathtubs, showers or water-filled reservoirs

This instruction manual can also be downloaded from our website at: [www.natureetdecouvertes.com](http://www.natureetdecouvertes.com)

## IMPORTANT INFORMATION

- Use of this appliance may cause discomfort or adverse effect to, persons suffering from skin disease or inflammation of the skin, persons with sensitive skin and persons having an allergic constitution. Stop using it if any discomfort occurs.
- This appliance is not to be used on the parts where the integrity of the skin is damaged (such as burns, blisters, lesions, incision, wounds, active skin disease or infection).
- This appliance shall not be used by, or on, persons wearing implants not compatible with its use.
- This appliance shall not be used by or on persons up to and including 14 years of age.
- **WARNING:** adverse effects may occur if the instructions are not followed.
- **WARNING:** prolonged or repeated application on one site may potentially damage the skin.
- **WARNING:** do not bypass or remove the inherent safety device of the appliance.
- Do not use on certain parts of the body such as eyes and the genital area.
- Only use disposable and consumable materials recommended by the manufacturer.
- Deactivate the output of the appliance before placing or removing the electrodes to or from the skin.
- Do not use the appliance in conjunction with flammable based cosmetic products, unless they are recommended by the manufacturer.
- Do not use the appliance near the eyes of the person on whom the appliance is used.
- Do not use the appliance on persons with an implanted electronic device (for example a cardiac pacemaker).
- Do not use the appliance close to the region of the heart.
- This device is designed for cosmetic use only and for individuals in good health. If this is not the case for you, please consult your doctor before use.
- **DO NOT** use over swollen, infected, or inflamed areas or skin eruptions, or in case of phlebitis, thrombophlebitis, broken capillaries, varicose veins, etc.
- **DO NOT** apply stimulation across your chest as the introduction of electrical current into the chest may cause rhythmic disturbances to your heart, which could be fatal.
- **DO NOT** use while driving, operating machinery, or during any activity in which involuntary muscle contractions may put the user at undue risk of injury.
- Never operate this device if it is not working properly, if it has been dropped or damaged, or submerged in water.
- Never drop or insert any object into any openings on the device.
- **DO NOT** use the device when it is charging.
- Use of controls or settings, or performing procedures other than those specified herein may result in injury.
- If the device begins to overheat, turn off immediately.
- Stop using the device immediately if you suspect it is malfunctioning, or if electric shock, fire or injury occurs.
- There are no user- serviceable parts inside the device.
- **Note:** disassembly will void the warranty.
- If using the device causes any discomfort, discontinue use immediately.
- If you suspect that you had a negative reaction from using the device, consult your doctor immediately.
- If you have suspected or diagnosed heart disease, you should follow the precautions recommended by your doctor.
- Caution should be used in the following situations:
  - for patients with suspected or diagnosed heart problems,
  - for patients with suspected or diagnosed epilepsy,
  - in cases of acute trauma or fractures with a tendency to hemorrhage,

- over areas of the skin which lack normal sensation.
- The long-term effects of electrostimulation are unknown. The safety of this device for use during pregnancy has not been established.
- Do not store the device in the sun or on a hot surface. High temperatures may be damaging. The device should be operated, transported, and stored in the manner specified in this manual. Failure to use and maintain the device in accordance with the instructions in this manual will void the warranty.
- Do not scratch or damage the device.
- You should stop using the device and should consult with your physician if you experience adverse reactions from the device.
- Do NOT immerse the device unit in water and do not put it in a dishwasher.
- DO NOT open or dismantle the housing of the unit, as there are no user serviceable parts inside.
- DO NOT drop or apply pressure to the outside of the unit since this may damage or weaken the bulb and other components.
- The appliance is not intended for use by children or persons whose physical or mental condition does not permit the safe use of the appliance.
- Do not place the device on hot heat sources (radiator, stove, fire, etc.).

## FUNCTION INTRODUCTION

The gua sha led for face and neck is an over-the-counter hand-held device designed to promote skin vitality and enhance the appearance of the face and neck. It combines 7 colors LED light, vibration and gentle heat to help enhance skin appearance and support a healthy-looking complexion.

### (A) Vibration

The vibration function has three intensity levels, each corresponding to a different massage mode.

- Level 1 corresponds to the CLEAN mode, which helps to remove impurities and purify the skin.
- Level 2 corresponds to the EMS mode, to gently stimulate the facial muscles.
- Level 3 corresponds to the HOT mode, to promote better penetration of skincare products using gentle heat.

### (B) Gentle heat

The gentle heat promotes the penetration of skincare products and makes their application more pleasant.

### (C) 7 Color LED Light

- GREEN LIGHT: supports skin's natural renewal process and helps maintain a youthful appearance.
- RED LIGHT: helps to improve the appearance of wrinkles and fine lines.
- BLUE LIGHT: helps soothe inflammation and calm sensitive skin.
- PURPLE LIGHT: helps promote lymphatic drainage.
- YELLOW LIGHT: helps reduce melanin production.
- CYAN LIGHT: helps reduce skin tension and the appearance of swollen capillaries.
- WHITE LIGHT: promotes skin relaxation and aids in the natural regeneration of skin cells.

## PRODUCT DRAWING

## PRECAUTIONS BEFORE USE

- If you wear contact lenses, remove them out before use.
- Wash your face to remove dirt such as makeup and sebum and apply lotion or serum to your skin after cleansing.
- Never use if there are any inflammation or abnormalities on your skin due to acne or pimples.
- Do not use it continuously on the same part of the skin but move it while using it. If it feels too hot, stop using it immediately.
- The vibration function remains active whenever the device is in contact with the skin.
- Do not use it on sensitive areas of the body such as the eyes or the genital area to avoid any potential injury or accidents.

## HOW TO USE

- Press the power button for 2 seconds to turn on the device.
- Press the button to select the function mode and the led color according to your needs. Each massage mode corresponds exclusively to specific color light functions:
  - Green, purple and white lights can only be used with CLEAN mode.
  - Red and yellow lights can only be used with EMS mode.
  - Blue and light blue lights can only be used with HOT mode.
- Place the massage head against your skin and start massage from bottom to top.
- Press the switch button for 2 seconds to turn off the device after using it.


### Note:

- The device can't be used during charging.
- The massage head must touch the user's skin for the vibration to start.
- The device will automatically turn off after 5 minutes.

## TIPS FOR USE

1. Massage from neck to chin.

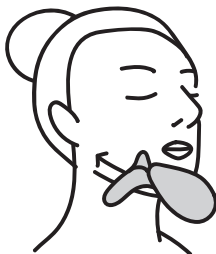2. Massage from chin to cheek.

3. Massage the forehead eyebrow to hairline.

It is recommended to use each mode for about 5 minutes, with 2 to 3 cycles per mode. Continuous use on the same area is not advised.

## BATTERY CHARGING

The device has a built-in rechargeable battery. The charging port is located at the bottom of device.

### Charging procedure:

- For charging, use the cable included in the box.
- Connect the USB-C inlet to the USB-A adapter output.
- The battery indicator flashes red during the charging process.
- Once the battery is charged, the red light will turn off.
- The maximum charging time is about 3 hours.
- When charging, please ensure that the appliance is placed within the safe temperature range of the battery (i.e. 0°C to 40°C).

## CLEANING AND MAINTENANCE

After treatment, always ensure the device is thoroughly cleaned. Never allow moisture to penetrate directly into the device. Ideally, use cotton or dry cloth for cleaning. Before using the device again, there should be no residue of gel, cosmetic cream, or serum on the device. Finally, always check that the device head is completely dry. Never use any sharp objects that could damage the surface of the applicator's head.

## TECHNICAL DATE

- Input: 5V  $\overline{\text{---}}$  1A
- Rechargeable lithium-ion battery
- Maximum charging time: approx. 3 hours
- Maximum operating time: approx. 1 hour
- Maximum EMS output (DC): 0.7mA/24V
- Dimensions: approx. H. 12 cm x W. 4.5 cm x L. 9 cm
- Weight: 95 g
- Battery capacity: 500 mAh

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