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PLEASE CAREFULLY READ THIS MANUAL AND KEEP IT IN A SAFE PLACE

WARNING

To mitigate its impacts and possible harm, do not use the product with the following medical electronic devices:

- Implanted medical electronic devices such as pacemakers.
- Life-sustaining medical electronic devices such as artificial hearts and lungs.
- Medical electronic devices such as electrocardiograms.
- Use is prohibited in flammable and explosive areas such as gas stations, high-dust workshops, and gas stations.
- Persons with neck ailments should consult healthcare professionals to avoid delaying treatment.
- This product is not a medical device and does not diagnose or treat neck disorders.
- Do not use this product while driving, climbing, or engaging in other dangerous activities, as it may lead to accidents.
- This device is intended for use on the neck area only.
- The heating function of this appliance should not be used for more than 15 minutes, as there is a risk of moderate heat burns. Do not use while asleep or in an unconscious state. Users who are sensitive to heat should use a lower heat setting and strictly control the use time.
- Individuals with circulatory disorders such as cerebral vascular diseases, inflamed or swollen neck skin, suppuration, or neck-related ailments should not use this product, as it may worsen symptoms.
- Keep this product out of reach of infants, children, or individuals who are entirely or partially unconscious. Do not let them use it, as it may lead to accidents.
- Individuals over 70 years old are prohibited from using this product. Please select a lower setting and strictly control the usage duration since those over 65 have more delicate skin and less sensitive nerve endings.
- Gradually increase the massage setting according to personal comfort; otherwise, the discomfort from a higher setting might cause muscle strain.
- This is a personal care massage product and cannot be used to substitute medicines or medical devices.
- Children and pregnant women are prohibited from using this product.
- Use this product under a doctor's supervision if you have sensitive skin, facial allergic dermatitis, or facial areas treated with cosmetic products.
- Users allergic to metal or heat must take precautions.
- This product has a heating function; please be cautious if you are insensitive to heat.
- If discomfort occurs during use, please stop using immediately.
- Do not disassemble this product yourself.
- Do not attempt to replace the battery yourself.
- If you have a neck injury, consult a healthcare professional before using this product.
- Use within the following operating environment: Temperature: 0°C to 40°C; Relative Humidity: <80%RH, no condensation; Atmospheric Pressure: 50Kpa to 106Kpa.
- Transportation and storage: Temperature Range: -20°C to +55°C; Relative Humidity Range: ≤95%RH, no condensation; Atmospheric Pressure: 50Kpa to 106Kpa.
- This product contains a lithium battery. Do not disassemble, modify, immerse in water, or expose to open flame, as it may result in damage or explosion. Do not place it in a high-temperature environment.
- Do not press on the heart when using, as this may

cause an accident or physical discomfort.

- Do not use this product with medical devices or with more than two of this product simultaneously, as it may cause nausea or physical discomfort.
- Do not use for purposes other than those specified in the instruction manual, as it may cause accidents.
- Do not disassemble, repair, or modify this product yourself, as it may lead to a malfunction or an electric shock accident.
- Do not use to massage parts of the body other than the neck.
- If a person simultaneously uses high-frequency surgical equipment and this massager, it may cause burns at the massager's electrode location and may also damage this massager.
- If this product is used near shortwave or microwave therapy equipment (within 1 meter), the output of this massager may be unstable.
- Electrode pads should avoid contact with wounds and scars during massage.
- While massaging, the electrode pads must not be placed near the heart area of the chest, as this could increase the risk of cardiac fibrillation.
- If you experience skin irritation when using the body massager, please refrain from using it, and wait for two to three days to see if your symptoms ease. If symptoms persist or worsen, please consult a physician.
- If you have eczema, allergies or asthma, you may be more likely to experience a skin irritation or allergy from using the massager.
- The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- Use of this appliance may cause discomfort or adverse effects to persons suffering from skin disease or inflammation of the skin, persons with sensitive skin and persons having an allergic constitution. Stop using it if any discomfort occurs.
- This appliance is not to be used on the parts where the integrity of the skin is damaged (such as burns, blisters, lesions, incision, wounds, active skin disease or infection).
- This appliance shall not be used by, or on, persons wearing implants not compatible with its use.
- This appliance shall not be used by or on persons up to and including 14 years of age.
- **WARNING:** adverse effects may occur if the instructions are not followed.
- **WARNING:** prolonged or repeated application on one site may potentially damage the skin.
- **WARNING:** do not bypass or remove the inherent safety device of the appliance.
- Deactivate the output of the appliance before placing or removing the electrodes to or from the skin.
- Do not use the appliance in conjunction with flammable based cosmetic products, unless they are recommended by the manufacturer.
- Do not use the appliance on persons with an implanted electronic device (for example a cardiac pacemaker).
- Do not use the appliance close to the region of the heart.
- The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- At end of life, the battery shall be properly disposed of.
- It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

In the following situations, please consult with a physician before use, otherwise accidents or discomfort may occur:

- Persons currently under medical treatment and experiencing abnormal physical symptoms.
- Persons with infectious diseases.
- Persons allergic to stainless steel electrode pads.
- Persons with heart or brain nerve abnormalities.
- Persons with sensory impairments or skin abnormalities.
- Persons with a body temperature above 38 °C (during a fever).
- Persons taking medications.
- Persons with abnormal blood pressure.
- Persons who are unable to express their own will freely.
- Persons with perception disorders due to severe peripheral circulation problems caused by conditions such as diabetes.

The device must not be used in the following situations:

- Persons with malignant tumors.
- Persons with a tendency to bleed.
- Pregnant women and women in childbirth.
- Users of cardiac pacemakers.
- Persons with acute suppurative inflammation.
- Individuals with severe cognitive impairments.
- Persons with metal implants.
- Persons with spastic paralysis.
- Individuals with clotting disorders.
- Individuals with abnormal nervous stress.
- Near the human heart.
- Individuals in a spasm phase following a cerebrovascular accident.
- Individuals with unclear consciousness, and limb and joint contracture deformities.

NOTE

- Please read the instruction manual thoroughly before using this product.
- If you feel physical discomfort due to abnormalities in this product, please stop using it immediately. Consult a doctor if you are experiencing unusual sensations in your body or skin, and follow medical advice.
- Do not use this device on top of other electronic devices, as it may cause malfunctions (for example, disrupting the timing of an electronic watch).
- If a fault occurs and the device cannot operate, immediately disconnect the power and remove it; otherwise, it may cause overheating or short circuits, resulting in an accident. Contact customer service for after-sales consultation.
- Do not plug or unplug the charging cable with wet hands, as it may cause electric shock or accidents.
- Do not use while sleeping, as it may cause the device to malfunction and the electrode pads to move to other parts of the body, causing discomfort.
- Please wear this product correctly on the neck, placing it symmetrically on both sides to avoid accidents or discomfort.
- When disconnecting the charging cable from the device, hold the plug rather than the cable, since pulling may cause an electric shock or damage the product.
- It is normal for the product to feel warm on the surface after prolonged use.
- When the product is working, muscles will contract without conscious effort: do not massage for an extended period of time. If discomfort is felt, take a break before using again.
- Do not forcefully fold the product when using or storing it to prevent damaging its components.
- When placing the product around the neck, take care not to pinch your hair or skin.

Please do not use or store this product in the following locations:

- In places exposed to direct sunlight
- In places with water

- In strong electromagnetic fields
- In hot and humid locations
- In places with excessive dust
- In locations exposed to vibration and shock
- Avoid placement near open flames.

- Keep away from hard objects, and avoid impacts and scratches.
- Clean the product with a slightly damp or soft cloth.
- Charge the product once a month if it is not in use for a long time and keep the battery fully charged when using the product.

Read operator's manual

WARNING: Do not use this appliance near bathtubs, showers or water-filled reservoirs.

RECOMMENDATION

- Prior to use, please ensure the main host buttons are working properly.
- If using the product for the first time or after the battery has been depleted from long-term storage, charge the product according to the instruction manual to ensure it works properly.
- Please dispose of this product in line with applicable city environmental protection standards to prevent environmental pollution.
- Recommended for use by individuals aged 16 and above.

PRECAUTIONS DURING HANDLING OF BATTERIES

- When used correctly, lithium batteries provide a safe and dependable source of power.
- However, if they are misused or handled improperly, leakage, venting, or in extreme cases, explosion and/or fire may occur.
- Keep batteries out of the reach of children.
- Do not force discharge batteries.
- Do not mix new and used batteries or batteries of different types or brands.
- Exhausted batteries should be immediately removed from equipment and properly disposed of.
- Do not heat batteries.
- Do not weld or solder directly to batteries.
- Do not dismantle batteries.
- Do not deform batteries.
- Do not dispose of batteries in fire.
- A lithium battery with a damaged container should not be exposed to water.
- Do not encapsulate and/or modify batteries.
- Store unused batteries in their original packaging away from metal objects. If already unpacked, do not mix or jumble batteries.
- Remove batteries from equipment if it is not to be used for an extended period of time unless it is for emergency purposes.

WARNING

KEEP OUT OF REACH OF CHILDREN


Swallowing can lead to chemical burns, perforation of soft tissue, and death. Severe burns can occur within 2 hours of ingestion. Seek medical attention immediately.

- Do not allow children to replace batteries without adult supervision.
- Always insert batteries correctly with regard to polarity (+ and -) marked on the battery and the equipment.
- Do not short-circuit batteries.
- Do not charge batteries.


DISPOSAL OF BATTERIES

- Batteries may be disposed of via communal refuse arrangements provided no local rules to the contrary exist.
- During transport, storage and handling for disposal, the following precautions should be considered:
 - Do not dismantle batteries.
 - Do not dispose of batteries in fire except under conditions of approved and controlled incineration.
- Store collected batteries in a clean and dry environment out of direct sunlight and away from extreme heat.
- Store collected batteries in a well-ventilated area.
- Do not mix collected batteries with other materials.
- Protect battery terminals.

PRODUCT DIAGRAM


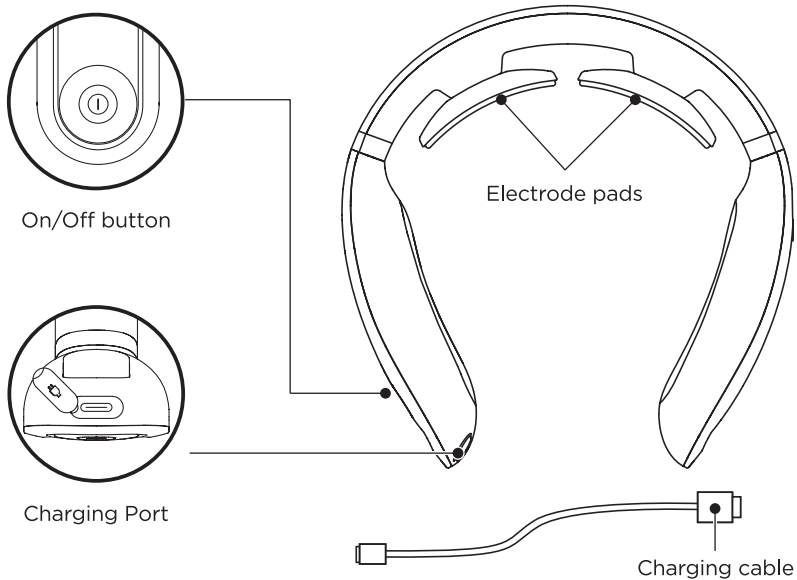
REMOTE CONTROL

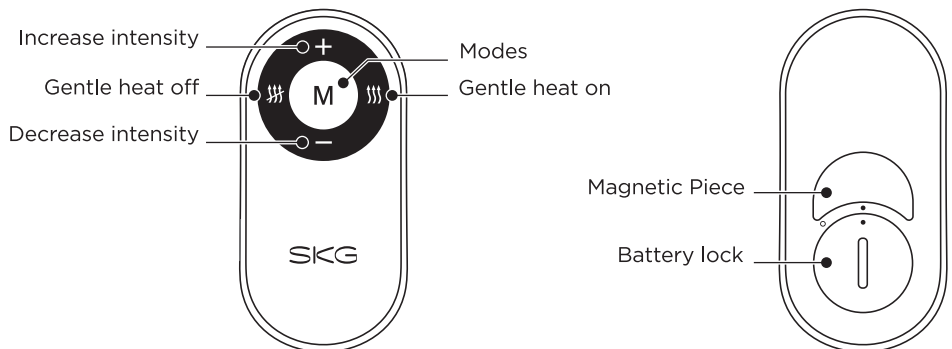PRODUCT PARAMETERS

Model	4098E	Rated Input Voltage	5V $\overline{=}$
Battery Capacity	1500 mAh	Rated Input Power	5 W
Transmitter frequency range	433,050 MHz – 434,790 MHz	Maximum radio-frequency power transmitted	20 mW
Power consumption in off mode	0,15W	The period after which the equipment reaches automatically off mode	<20 minutes

TIPS FOR USE

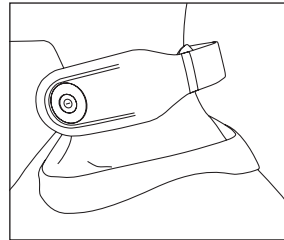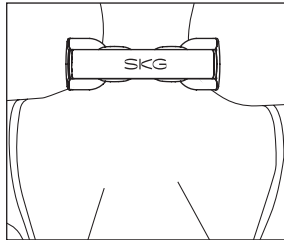1. Moisten the skin on the cervical spine area with a damp towel.
2. Wear the product around the neck, tightly against the skin.
3. Ensure the electrode pads are firmly in contact with the skin.
4. After adjusting to the proper position, turn on the device and set the massage intensity and mode to a comfortable level.

DEVICE OPERATION

Power On

- Press and hold the «①» button for 2 seconds (the buzzer beeps once, and the indicator light turns on), entering the massage mode and intensity level 1.

Intensity Adjustment

- When the device is on, press the «①» button once each time: the buzzer beeps once, and the intensity increases by one level (levels 1-10).

- When the intensity reaches level 10, pressing the button again will cause the buzzer to beep twice, and the intensity will not increase further.

Power Off

- Long press the «①» button for 2 seconds to power off (the buzzer will beep twice).

REMOTE CONTROL OPERATION

Mode switching

- Press the «M» button to switch modes: the buzzer will beep once for massage mode 1, twice for massage mode 2, and three times for massage mode 3.
- When mode changes, intensity will be reseted to level 1.
- One massage cycle lasts 10 minutes and then ends automatically.

Intensity adjustment

- «+» button: increases intensity
 - Press once: a beep sounds and the intensity increases by one level (levels 1 to 10).
 - When the intensity reaches level 10, pressing '+' again triggers two beeps; the intensity remains at level 10.

- «-» button: decreases intensity
 - Press once: a beep sounds and the intensity decreases by one level (levels 10 to 1).
 - When the intensity reaches level 1, pressing '-' again triggers two beeps; the intensity remains at level 1.

The gentle heat function

- Press the «|||» button to turn on the gentle heat function: the buzzer beeps once, and the electrode pads of the device heat up, with the temperature stabilizing between 37-43 °C.
- Press the «|||» button to turn off the gentle heat function: the buzzer beeps twice, and the heating function stops.

MODE EXPLANATIONS

Mode 1

- Longer electrical pulse duration ~5 seconds (5 seconds of massage -> rest -> 5 seconds of massage, etc.)
- Provides gentle, continuous pulses to help relax the muscles.
- Suitable for daily relaxation and mild neck discomfort, relieving a sense of fatigue.

Mode 2

- Shorter electrical pulse duration ~1 second (1 second of massage -> rest -> 1 second of massage, etc.)

- Especially suitable for relieving neck tension caused by long hours of work or study.
- Delivers stronger muscle stimulation to effectively relax stiff muscles.

Mode 3

- Reproduces a tapping effect, similar to the techniques used by professional masseurs.
- Helps stimulate blood circulation.
- Provides dynamic and comprehensive stimulation of the cervical region, particularly suitable in cases of intense fatigue.

TROUBLESHOOTING INSTRUCTIONS

If the product is not working properly, please check for the following causes.

If our suggestions do not solve the problem, please return the product to the dealer for repair.

Condition	Possible Cause	Solution
Unable to start	Low battery power	Charge the product
	The product is not turned on	Long-press the power button to turn on
	The product is still connected to the household power socket	Remove the charging cable from the device
Unable to charge or no charging light displayed	Charging cable is not properly inserted into the power socket	Replug the charging cable
	Wrong charging cable used	Use a Type-C connector cable

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S.A. Au capital de 37 653 500 €