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FR

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AIR COMPRESSION FOOT MASSAGE**EN**

Warning: Please read this instruction manual carefully before using the instrument and keep it for future use.

A. GENERAL SAFETY INSTRUCTIONS

1. The device is intended exclusively for private and not for commercial use. Only use this device as describe in this instruction manual.
2. Before to use, make sure that the supply voltage corresponds with the voltage marked on the rating label.
3. Do not leave the device unattended when plugged in.
4. Do not handle the unit or plug with wet hands.
5. Unplug from the power outlet after use and allow the appliance to cool down completely before to clean it or to move it.
6. When disconnecting the cable, always pull at the plug; do not pull on the cable itself.
7. When using the device, always pull the power cord completely from its housing to prevent overheating.
8. If the cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid an electric shock.
9. Never try to repair it by yourself. It must be repaired by the manufacturer or its service agent or a similarly qualified person in order to avoid an electric shock.
10. The power cable must be checked regularly to avoid damage. If the power cable is damaged it must be replaced by a qualified service center to avoid any danger.
11. Always place your device on a dry surface.
12. Do not use outdoors.
13. Do not let the power cord on the edge of a table or work surface, and prevent it comes into contact with hot surfaces.
14. Keep it away from heat sources such as radiators, to avoid deformation of plastic parts.
15. Do not touch the hot surface. Keep the appliance out of reach of children.
16. This unit is not designed to be operated by means of an external timer or a separate remote control system.

B. SPECIAL SAFETY INSTRUCTIONS

1. This device must not be used by children under 3 years of age as they are unable to react in the event of it overheating.
2. This device must not be used by small children (3-8 years) unless a parent or guardian has pre-set the necessary controls, or the child has been given instructions on how to use the heating pad safely.
3. This device can be used by children aged 8 or over, and by persons with physical, sensory or mental disability or with limited experience or knowledge, provided they are under supervision or have received instructions on how to use it.
4. Children must not play with it.
5. Children must not be allowed to clean or maintain this device without supervision.
6. Make sure that you do not fall asleep while using this device!
7. If you are very tired while using it, you can select level 1 for your safety.
8. You must examine the device frequently for signs of wear and tear or damage. If you detect signs of wear or tear or damage, or if the device has been used inappropriately stop to use it.
9. Do not attach safety pins, needles or pins of any kind to this device.
10. This product is not intended for use in hospitals, clinics etc
11. This is not a medical product.
12. Do not use with a helpless person, child or any person with low sensitivity to heat.

C. PRECAUTIONS


1. People with the following conditions should consult their physician before using the product:
 - Who uses a medical electronic device, which is embedded in the body, such as a pacemaker.
 - Who suffers from heart disease.
 - Who suffers from backbone injury or disease.
 - Who is suffering from bone fragility, especially on the backbone.
 - After surgery.
 - Who is pregnant or in a period of monthly illness.
 - Who suffers from malignant tumors.
 - Who is suffering from osteoporosis.
 - Who is suffering from blood circulation problems or thrombosis.
 - Who is suffering from skin disease or infection.
 - Who is undergoing medical treatment or feels some physical abnormality previously or is being treated now.
 - Who has some abnormality or deformity of limbs or backbone.
 - Who has injuries, diseases due to illness or accident.

D. INTRODUCTION

1. Leg massager
2. Air pressure hose
3. LCD screen
4. Control panel

E. HOW TO USE IT

- Slip on each massager, keeping your feet flat on the bottom. Note that each is marked L and R for each foot.
- Use the strap to tightly secure your foot into the massager
- Use the strap to tightly secure your leg into the massager
- Move your foot to make sure that the massage is securely fastened
- Insert the two air pressure connections into the controller
- Press the POWER button ❶ and the LED will light on.
- Choose the air pressure strength by pressing POWER button shortly : low – normal – high intensity
- Press MODE button ❷ to choose from the 3 different massage modes:
 - I. Calf + foot
 - II. Calf only
 - III. Foot only
- Press HEAT button ❸ to start heating (on foot only) and select low – normal – high heat intensity
- The device will turn off automatically after 20 minutes
- If you want to turn it off before, press POWER button ❶ for 3 seconds.

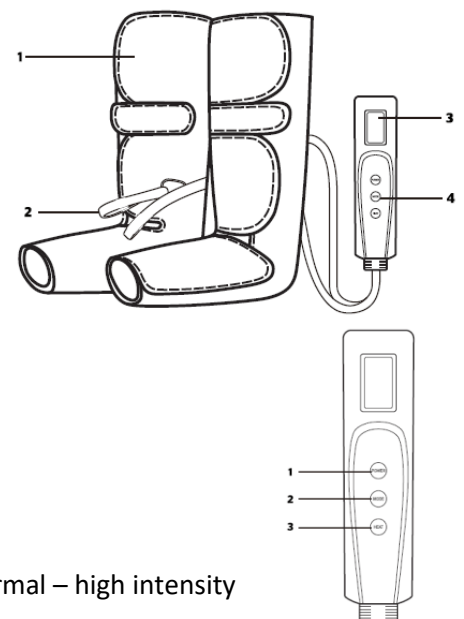F. CARE AND MAINTENANCE

Clean the product only with a dry or slightly damp, lint-free cloth.

Do not use abrasive cleaners to clean the unit.

Store in a cool, dry place, out of direct sunlight. Avoid humidity.

G. TROUBLESHOOTING

Problem	Cause	Solution
The device does not work.	The unit is not connected to a power source.	Make sure the unit is fully connected to the power.
	The power button was not pressed.	Press the power button.
The device suddenly stopped working.	The unit was running more than 15 minutes.	Restart the unit to resume use.

H. FEATURES

- Material: polyester, ABS
- Power: AC 110V-220V DC 12V 1.5A (adaptor 180 cm included)
- Massage : 3 modes (foot and calf / calf / foot)
- 3 intensity (25 KPa / 35-40 KPa / 45-48 KPa)
- Heating function (on foot only) : ON/OFF + 3 heating levels

The European directive 2012/19/EC about electronic and electric waste, requires that you can't throw away defective domestic appliance with common waste. Used device has to be collect separately to optimize the recuperation rate and the recycling of the materials to reduce the impact on the health and the environment.

Importé par / Importiert von / Geïmporteerd door / Imported by :
CAPTELEC, 59170 Croix - FRANCE