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Warning: Please read this instruction manual carefully before using the instrument and keep it for future use.

A. SAFETY INSTRUCTIONS

1. This appliance is intended for household use only. Use the appliance only as described in this manual.
2. Before plugging in the appliance, make sure that the electrical voltage in your home corresponds to that indicated on the appliance's rating plate.
3. To prevent electric shock, never immerse the power cord, plug or any other part of the appliance in water or any other liquid.
4. Do not handle the appliance or the plug with wet hands.
5. When unplugging the cable, always pull it by the plug, do not pull on the cable itself.
6. When using the appliance, always unwind the power cord completely from its housing to avoid overheating the appliance.
7. Never attempt to repair the appliance yourself. Always have it repaired by a service center or a qualified repair center.
8. Do not let the power cord hang over the edge of a table or work surface, or allow it to come into contact with hot surfaces.
9. Keep the device away from heat sources such as radiators, to avoid deformation of the plastic parts.
10. Do not use abrasive cleaning products to clean the appliance.

B. SPECIFIC SAFETY INSTRUCTION


1. Do not let water enter the product or be splashed near the surface of the product.
2. Do not operate the button with wet hands.
3. Keep the buttons from water or other liquid.
4. Water contact with electricity is dangerous. Do not leave the product in a humid environment (such as in the bathroom or shower room or swimming pool).
5. Before connecting the product to the power supply, check that the product voltage and supply voltage are the same.
6. If the plug, power cord or product itself is damaged, do not use the product.
7. In order to prevent the risk of electric shock, do not open the product outer cover.
8. Product is not intended for use by the physically or mentally impaired (including children), unless there is supervision or guidance by the people responsible for their safety when they use the product.
9. Children should not use this product. They should be under supervision to ensure that they will not treat this product as a toy.
10. If an unusual sound or noise is being heard from the product, turn off the product and immediately disconnect the power connection.
11. Do not use this product in high temperature conditions or beside hot objects such as heaters, stoves.
12. Do not apply medical or cosmetic creams on the product or on your skin before using this product.
13. Do not destroy, bend stretch, twist or cut the power cord to avoid electric shock.
14. People who are sensitive to heat should not use the heat function of product.
15. Please consult your doctor before using if you have the following situations
 - Receiving medical treatment or experiencing physical discomfort
 - Suffering from hypertension, heart disease, cerebral arteriosclerosis and cerebral thrombosis
 - Suffering from cancer
 - Pregnant or in a period
 - Suffering from osteoporosis
 - Using a pacemaker or implantable medical devices, artificial joints

- Suffering from acute diseases
- Minors

16. Do not fall asleep when using the product.
17. Do not use this product when you are drunk or feeling uncomfortable.
18. Do not directly contact the massager pad with your hands, feet joints, knee, chest and abdomen.
19. Do not use if the product's cloth cover has been damaged.
20. Do not use the product before going to bed.
21. Do not enable the heating function when the product is being used in a low room temperature. It is recommended that you should gradually adjust the temperature to abnormal range before enabling heating.
22. If the product is used immediately after being transferred from a cold environment to a warm environment, the accumulation of water vapor inside the product may affect the normal operation of the product or even cause mechanical failure. You are suggested to put the product in a normal room temperature for 1 hour before starting to use it after the transfer.

C. INTRODUCTION

- 1) Loading port
- 2) Indicator light
- 3) Intensity icon
- 4) Temperature icon
- 5) Vibration icon
- 6) Temperature selection icon
- 7) ON/OFF button and screen lock

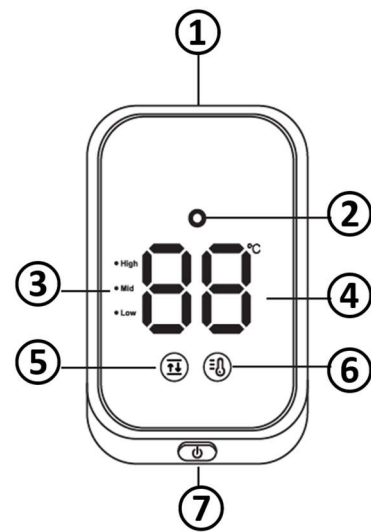D. HOW TO USE

1) Massage areas

ELBOW

KNEE

SHOULDERS

2) Function

- 1- Press the ON/OFF button (7) for 2 seconds to turn on the massager. A beep sound will be emitted and the indicator light will light up to indicate that the massager is on. If you do not perform any action for 10 seconds, 2 beeps will be emitted and the massager will turn off.
After 30 minutes of use, the massager will automatically turn off

- 2- Attach the massager using the straps to your elbow, knee or shoulders

- 3- Tap the vibration icon (5) to choose the intensity.

Press 1 time for low intensity, 2 times for medium intensity, 3 times for high intensity, 4 times to stop vibration (3).

- 4- Press the temperature selection icon (6) to choose the temperature.

Press the icon again to choose the desired temperature (4): 45°C; 50°C; 55°C; 60°C; 65°C

Press a 6th time to stop the heat.

5- Short press the ON/OFF button ⑦ to lock the screen.

The green light turning red indicates that the screen is locked.

Red: locked

Green: unlocked

Press the ON/OFF button ⑦ again to unlock the screen

3) COMMONLY ENCOUNTERED PROBLEMS

ISSUE	REASON	SOLUTION
Does not turn on	The product is plugged into a power outlet	Remove the charging cable from the product
	The battery is low	Charge the battery
	The ON/OFF button is not pressed firmly	Press the ON/OFF button firmly until it clicks into place.
The product does not heat up	The product is not for use outdoors, only at room temperature.	Use the product at a temperature between 15-35 degrees Celsius
The product stops during use	The battery is discharged	Charge the battery
	After 30 minutes of use, the product automatically shuts down	No problem with the product

4) MAINTENANCE ADVICE

Clean the product only with a dry or slightly damp, lint-free cloth.

Do not use abrasive cleaning products to clean the appliance.

5) TECHNICAL CHARACTERISTICS

- Materials: Synthetic textile, ABS
- Power supply: 1100mAh Li-Po battery, 5V, 7W USB charging cable included (power adapter not included)
- Massage: 3 intensity levels - Heat: 5 intensities (45-50-55-60 65°C)
- Autonomy: up to 6 hours depending on mode
- Charging time: 4-5 h

The European directive 2012/19/EC about electronic and electric waste, requires that you can't throw away defective domestic appliance with common waste. Used device has to be collect separately to optimize the recuperation rate and the recycling of the materials to reduce the impact on the health and the environment.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capacities, or by persons lacking in experience or knowledge, unless they are supervised or instructed on the use of this apparatus by a person responsible for their safety. Children should be supervised to ensure they do not play with this appliance.

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