LANAFORM

COMFORT

Heating Pillow


CORDLESS HEATED CUSHION

EN ENGLISH-04 CORDLESS HEATED CUSHION
FR FRANÇAIS-10 COUSSIN CHAUFFANT SANS FIL

NL NEDERLANDS-17 DRAADLOOS VERWARMD KUSSEN
DE DEUTSCH-24 KABELLOSES HEIZKISSEN



1

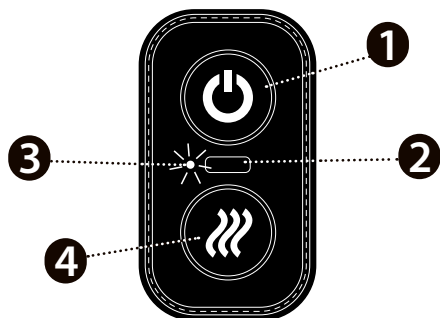2

3

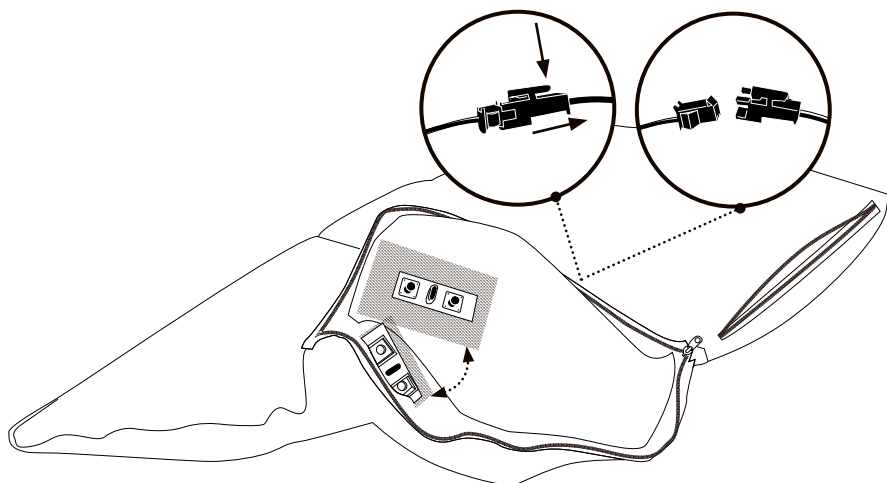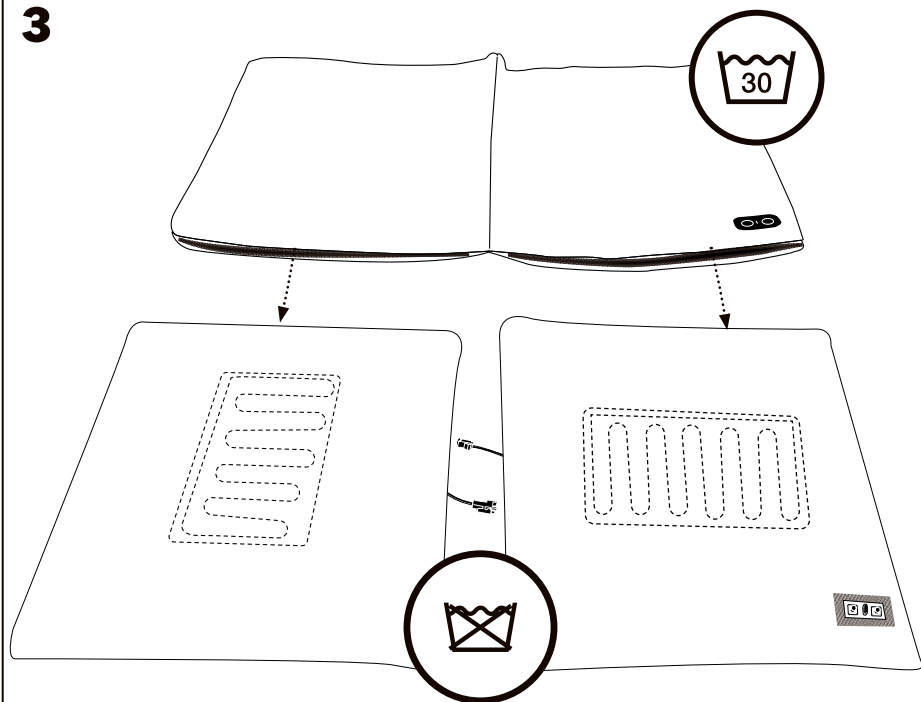Contents of the packaging

Check that the outer carton of the device is intact and all items are present. Before use, ensure that the device and accessories do not show any visible signs of damage and that the entire packaging has been removed. If in doubt, do not use the device and contact your dealer or the customer services department indicated.

- 1 Heating Pillow
- 1 USB cable
- 1 instruction manual

Explanation of symbols

The following symbols are used on the instructions for use, packaging and device nameplate:

	WARNING: Warning instruction indicating a risk of injury or danger to your health
	CAUTION: Safety note indicating possible damage to the device/accessories
	NOTE: Note on important information
	Dispose of in accordance with the WEEE (Waste Electrical and Electronic Equipment) European Directive
	Read the instructions
	The device is double insulated and therefore complies with Protection Class 2
	Do not pierce with needles
	Do not use when folded or rolled up
	Not for use by young children (0-3 years)
	This product meets the requirements of the applicable European and national directives.
	Maximum washing temperature 30°C, ultra-delicate cycle
	Do not bleach
	Do not tumble dry
	Do not iron

	Do not dry clean
	Manufacturer
	Separate the product and packaging components, and dispose of in accordance with municipal requirements.
	Separate the packaging components, and dispose of in accordance with municipal requirements.
	Paper packaging is either recyclable or made from recycled materials

Important instructions. Keep for later use

 **Please read all the instructions before using your product, in particular these few fundamental safety instructions.**

 **The photographs and other product images in these instructions and on the packaging are as accurate as possible but do not claim to be a perfect reproduction of the product.**

Warning

- **Only** use this device in accordance with the instructions given in this manual.
- **Before** plugging in your device, check that your mains voltage is the same as that of the device.
- **Children** under three years of age must not use this device due to their inability to react in the event of overheating.
- The device can be used by children between 3 and 8 years of age under supervision, with the control always set to the minimum temperature value.
- This device can be used by children aged 8 years and over and by people with reduced physical, sensory or mental capacities, or lack of experience and knowledge, providing they are supervised or have been instructed

on the safe use of the device and understand the risks involved. Children must not be allowed to play with the device. Cleaning and maintenance by the user must not be carried out by children without supervision.

- This heating cushion is not intended to be used in a hospital. This unit is intended only for family use.
- **Never** use this device in a room in which aerosols (sprays) are used or in a room in which oxygen is administered.
- **Do not** use the heating cushion with other heating devices or to heat other things (foot warmer, hot water bottle, etc.)
- If using the electric heating cushion for several hours, it is recommended to set the switch to the lowest temperature level to avoid overheating.
- **Do not** pinch the cushion (especially on beds with variable positions or on mattresses with separate elements).
- **Do not** use the heating cushion on animals.
- **Always** unplug the device after use or before cleaning it.

- To disconnect the device, switch off the device by pressing the On/Off button, then remove the USB cable if necessary, then remove the plug from the power socket.

- Disconnect the heating cushion if it becomes abnormally hot.

- If a device falls into water unplug it immediately, before retrieving it.

- When storing your heating cushion for a long time, store it in a dry place. Carefully check the condition of your heating cushion before operating it again for the first time. If the USB cable is damaged, it must be replaced with the same type of cable which is available from the supplier or its after-sales department.

- Do not use attachments that are not recommended by LANAFORM or supplied with this unit.

- Use your heating cushion with the included USB cable or equivalent.

- Do not carry this device by its USB cable or use the cable as a handle.

- An electrical device must never be left plugged in and unattended. Unplug it whenever you are not using it.

- Do not use this device if the plug is damaged, if the device is not working properly, if it has been dropped on the floor, if it is damaged or if it has fallen into water. If one of these things happen, have the device examined and repaired by the supplier or their after-sales department.

- Do not pull, twist, or bend the cables too strongly.

- Do not scratch the heating cushion

with a sharp object. Do not attach the heating cushion using drawing pins or other metal objects.

- Do not expose your heating cushion to sunlight and do not iron it. This will prevent you from damaging the insulation on the electrical circuit.

- Do not use the heating cushion when it is wet and do not turn it on to dry.

- Take every precaution to avoid burns.

- If you are concerned about your health, consult a doctor before using this device.

- If you experience any pain while using this device, stop using it immediately and consult your doctor.

- Never use this device on any part of the body that is swollen or inflamed or if you have a rash.

- When using the heating cushion, if you notice signs of abnormal redness on the skin, or even a feeling of discomfort or pain, stop using it for at least one hour. If the phenomenon persists, postpone use of the cushion and consult your doctor immediately, as necessary.

- Electrical and magnetic fields of this electrical device may interfere with the operation of pacemakers. However, they are well below the limits: electric field strength : Max. 5,000 V/m, magnetic field strength : Max. 80 A/m, magnetic induction : 0.1 mT max. Therefore, please ask your doctor and your pacemaker manufacturer for advice before using this electric heated blanket.

- The user may experience skin burns or thermal shock when falling

asleep with the device set to a high temperature.

- **When storing the device, allow to cool before folding.**
- **Inspect the device frequently for signs of wear or damage. If signs of wear or damage are found, or if the device has been used inappropriately or is not working, return to the supplier before switching on again.**

The device's cable may cause entanglement, choking, tripping or crushing hazards if not properly arranged. The user must ensure that the protruding connections and cable are arranged in a safe manner.

Presentation

Thank you for choosing the Heating Pillow by LANAFORM. Cozy and comfortable, this 2-in-1 heating cushion will delight even the chilliest of people! Thanks to its ability to heat quickly, the Heating Pillow gives you soothing warmth in just a few minutes, covering you in comfort without parallel. Equipped with two heating intensities and a rechargeable battery, the Heating Pillow lets you relax comfortably while staying warm, whether on your sofa as a heating cushion or as a heated seat cover.

Description 1

- 01 On/Off button
- 02 USB connector
- 03 On/Off & temperature indicator LED
- 04 Temperature setting

Preparation before use 1 & 2

The LANAFORM heating cushion is intended for warming the human body. It should be

used either:

As a cushion:

In this case, make sure that the cushion is closed using the zipper on the side of the product.

As a heated seat cover:

Unfold the cushion by opening the zipper around the rim of the cushion. The heating zones are located on the outer surfaces of the pillow (where the **On/Off** button is placed). When using the pillow as a seat cover, make sure to position it so that the **On/Off** button is against the part of your body that needs to be warmed.

Charging 2

You can charge your device using the provided USB cable. Connect the USB charging cable to a suitable power adapter and the device. Then plug the power adapter into an appropriate **power outlet (2.0 A)**.

While charging, the red indicator light pulses, meaning it slowly turns on and off repeatedly. The indicator light will turn green once the device is fully charged. Be sure not to turn on the pillow when you wish to charge it, as the charging process will not be recognized if the product is on.

Safety

This heating cushion is equipped with a safety system. The electronic sensor prevents the cushion from overheating over its entire surface by automatically shutting it off in case of failure.

The heating cushion is equipped with a safety over-temperature protection and is set to run continuously for 45 minutes.

Instructions

1. To turn on the device, press the **On/Off** button for 2 seconds. The green LED light comes on.
2. The heating cushion has 2 temperature intensity levels. Select the one you want (green LED = low intensity and red LED = high intensity) by successive-ly pressing the corresponding key.
3. Once the battery is charged, the heating cushion can run for approximately 90 minutes. It is also possible to use the heating cushion with an adapter. To do this, insert the USB cable into the connector and connect the adapter to the power outlet.

Warning

If the red or green indicator light flashes more quickly, it means the battery needs to be recharged! In this case, **turn off the device** by pressing the ON/OFF button and follow the charging steps outlined above.

Note: If the device is not turned off beforehand, it will continue flashing red/green and will not charge. If you wish to use the pillow while it is plugged into the mains, and the light is flashing quickly, the battery must be charged for 15 minutes before heating can begin.

If you are using the Heating Pillow using an adapter, please note that the USB port may become very hot. Please do not touch it while using the product

Cleaning and maintenance 3

When you are not using the heating cushion and before cleaning it, follow these instructions:

- Check that the appliance is turned off.
- If necessary, disconnect the adapter and the USB cable.
- Leave the heating cushion to cool before cleaning.

The cushion cover can be removed and machine washed at 30°C on a programme for delicate textiles.

- 01 Open the 2 zippers allowing access to the heating pads.
- 02 Then remove the Velcro from the control unit.
- 03 Disconnect the connection cable between the 2 heating pads.
- 04 Then remove the 2 heating pads. You can now machine wash the cover.
- 05 Once the cover has been machine washed and air dried, insert the 2 heating pads and connect them again.
- 06 Place the control unit securely in place, making sure the USB connector is correctly positioned.
- 07 Replace the zippers. You can now use the heating cushion

The cushion can be cleaned by hand with a damp cloth or sponge and water, or possibly with a little liquid detergent for delicate fabrics.

 **Be careful not to put the heating cushion in the machine or dry clean it. It should not be twisted, dried in a machine, calendered or ironed.**

When drying the cushion, do not hang it up with clothes pegs or similar. Do not expose the heating cushion to heat sources such as the sun, a radiator, etc.

 **Do not under any circumstances switch on the heating cushion to dry it.**

Storage

If you are not using the cushion for a long time,

it is advisable to store it in its original packaging in a dry place and not to place any objects on top.

Troubleshooting

⚠ Problem	🔍 Reason	😊 Solution
The cushion does not heat up	Battery discharged	Connect the USB cable to the cushion connector and to an adapter to recharge the battery or use the cushion by connecting the USB cable directly to the power socket.

Technical data

- Product reference: LA180504
- Voltage : DC5V 2.0A
- Power: 10 W
- Dimensions: 42 x 42 x 11cm

Waste disposal advice

The packaging consists entirely of environmentally safe materials that can be deposited at your local sorting centre to be used as secondary materials. Cardboard can be placed in a paper collection container. Packaging films should be dropped off at your local sorting and recycling centre.

If you are no longer using the device, dispose of in an environmentally friendly manner and in accordance with legal guidelines.

Limited warranty

LANAFORM guarantees that this product will be free from any material and workmanship defects for a period of two years from the date it was purchased, other than as specified below.

The LANAFORM warranty does not cover damage caused by normal product wear and tear. In addition, the warranty for this LANAFORM product does not cover any damage caused as a result of misuse or improper use, accidents, attaching any unauthorised accessories, modification of the product or any other condition of any kind beyond the control of LANAFORM.

LANAFORM shall not be held liable for any type of incidental, consequential or special damage.

All implied warranties of fitness for the product are limited to a period of two years from the original date of purchase, provided that a copy of the proof of purchase can be presented.

Upon receipt, LANAFORM will repair or replace your device, as appropriate, and return it to you. A warranty claim can only be made through the LANAFORM Service Centre. Any maintenance of this product entrusted to any person other than the LANAFORM Service Centre cancels this warranty.

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